


HQTS	INSPECTION REPORT Pre-Shipment Inspection	Report No: H-20062506
-------------	---	-----------------------

INSPECTION INFORMATION			
Client	Bunzl		
Supplier	Wuhan Careland Protection Co., Ltd		
Factory	Wuhan Huashida Protective Products Co., Ltd		
Inspection Instruction	GWI-H-004-HQT,SWI-1833-HQT-BUNZ		
P.O. No.	126759	5.1.1c	
Item No.	4.05.053.00	Sampling Plan	ISO2859 : Single,Normal
Inspection Date	09/06/2020	Inspection Level	G-II
Inspection Sequence	Original inspection	Sample Size	500 Pieces
Inspection Location	WuHan, China.	AQL	Critical:0, Major:2.5, Minor:4.0
Product Description	FACE MASK		

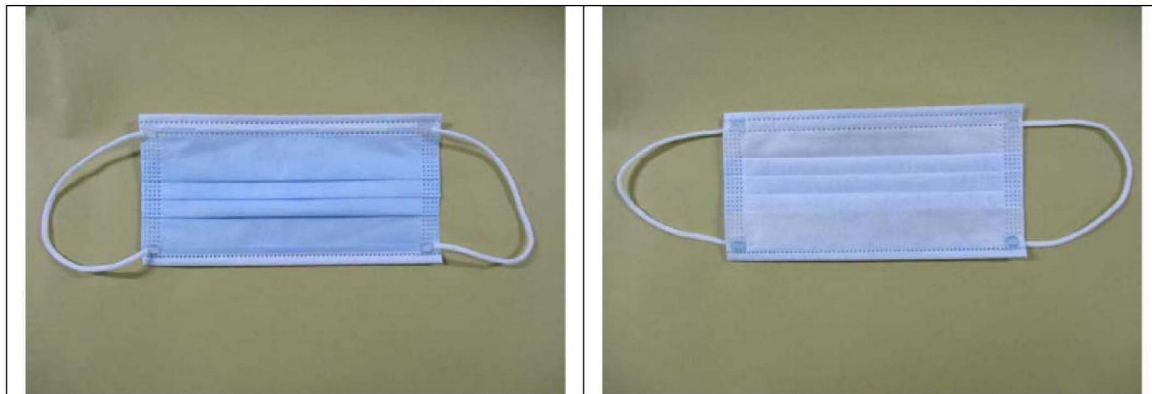

INSPECTION RESULT SUMMARY			
Category	Conform	Not Conform	Pending
1.Quantity	X		
2.Workmanship	X		
3.Style, Color & Documentation	X		
4.Field Test & Measurement			X
5.Shipping Mark & Packaging			X
OVERALL CONCLUSION	PENDING for client's evaluation		

Approved by HQTS Supervisor:	5.1.2e
-------------------------------------	--------

This report reflects the facts as recorded by HQTS at the time and place of inspection. It does not relieve the manufacturers from their contractual obligations nor prejudice client's right for compensation for any apparent and/or hidden defects not detected during our random inspection or occurring thereafter.

SPECIAL ATTENTION POINT (SAP)

1	Quantity	NIL
2	Workmanship	NIL
3	Style, Color & Documentation	NIL
4	Field Test & Measurement	
4.1	For microbial cleanliness test, no suitable apparatus was provided for checking by factory, but provided a test report for reference.	
5	Shipping Mark & Packaging	
5.1	The marking on master carton was different from spec, as photos.	
5.2	The certificate of qualification inside of color box was in Chinese.	



SAP 5.1 on spec



SAP 5.1 actually found



SAP 5.2 actually found

NA

INSPECTION DETAILS




5.1.1c


2. Workmanship	Conform		
	Critical	Major	Minor
Defects Description			
Poor welding of ear loop			4
Wrinkle mark on surface			3
Pilling on corner			1
Total Found:	0	0	8
Allowed:	0	21	21
Sample Size:	500	500	315

The following field test included in workmanship check

Labeling check

Defect Photos

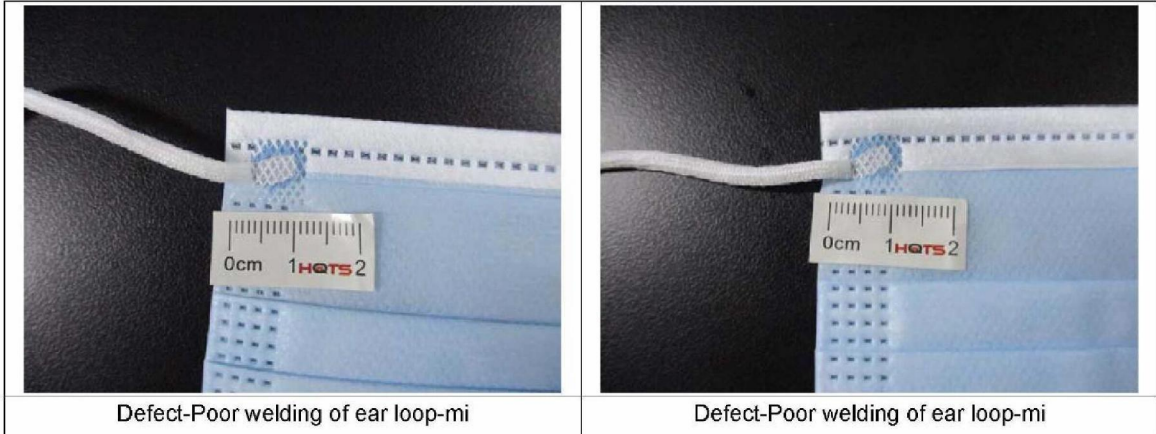

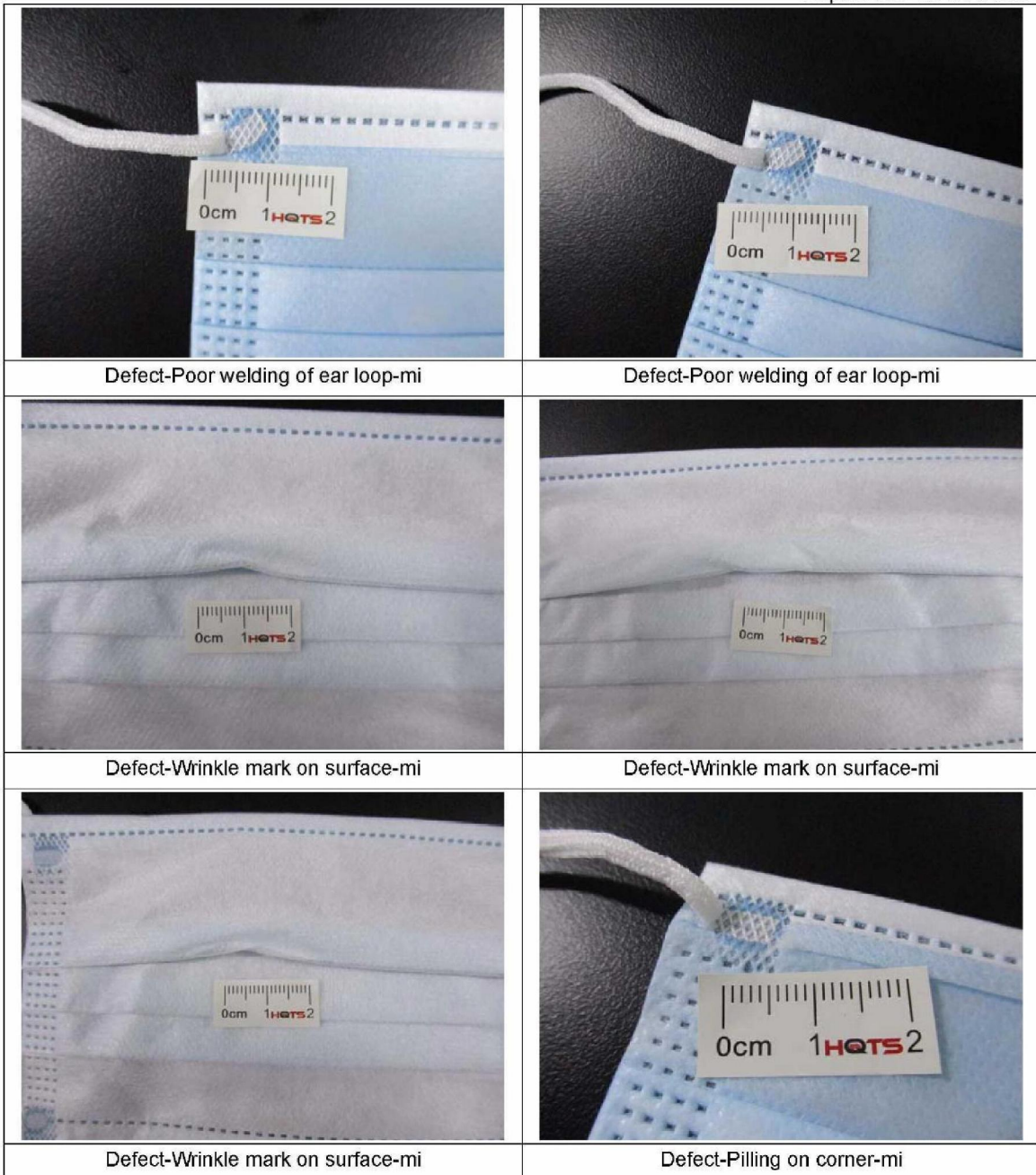

3. Style, Color & Documentation		Conform		
DESCRIPTION	Conform	Not Conform	Pending	
The inspected samples were consistent with	-	-	-	
Purchase Order	X			
Product Specification	X			
Instruction manual and Document Verification	-	-	-	

Report No: H-20062506

Client Comment	-	-	-
-----------------------	---	---	---

4. Field Test & Measurement				Pending		
Item	Specification	Sample Size	Conform	Not Conform	Pending	
4.1	Carton gross weight check	1 carton	X			
4.2	Smell check	5 samples	X			
4.3	Carton dimension check	1 carton	X			
4.4	Labeling check	5 samples	X			
4.5	Quantity per carton check	All selected cartons	X			
4.6	Barcode Readability Check	5 samples	X			
4.7	Product dimension check	5 samples	X			
4.8	Product weight check	5 samples	X			
4.9	Fabric weight check	5 samples/layer	X			
4.10	Transportation drop check	5 samples	X			
4.11	Fitting test	5 samples	X			
4.12	Tension test of ear loop	Min.10N	5 samples	X		
4.13	Bacterial filtration efficiency(BFE) test	BFE \geq 98%,Delta P <60	5 samples	X		
4.14	Differential pressure test	<60	5 samples	X		
4.15	Synthetic blood penetration test	\geq 16	5 samples	X		
4.16	Microbial cleanliness test				X	

Checking Item	Specification	Result
Fabric weight check	External layer: 28 gsm (+/-2gsm) Middle layer: 25 gsm (+/-2gsm) Internal layer: 23 gsm (+/-2gsm)	External layer: 27;27;27;27;27; Middle layer: 24;24;25;24;24; Internal layer: 22;23;21;23;22;

Item No.	Carton dimension (cm +/-5%)		Carton gross weight (kg +/-10%)	
	Specification	Result	Specification	Result
4.05.053.00	52x39x38	52x40x39	8.8	8.15

Item No.	Color box dimension (cm +/-5%)		Color box gross weight (g +/-10%)	
	Specification	Result	Specification	Result
4.05.053.00	-	19x10x9; 19x10x9; 19x10x9; 19x10x9; 19x10x9;	-	181.3; 183.59; 183.01; 183.65; 183.27;

Item No.	Product dimension (cm +/-5mm)		Product weight (kg +/-10%)	
	Specification	Result	Specification	Result
4.05.053.00	17.5x9	17.5x9;	-	2.93;

Report No: H-20062506


		17.5x9; 17.5x9; 17.5x9; 17.5x9;		2.90; 2.92; 2.92; 2.91;
--	--	--	--	----------------------------------

5. Shipping Mark & Packaging	Pending
---	----------------

5.1 Shipping mark photos



5.2 Packaging photos

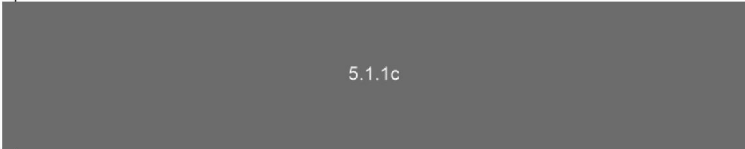







6. Reference Note	
6.1	Bacterial Filtration Efficiency (BFE) Test, the result was automatically converted result by the Particle filtration efficiency (PFE) , since no corresponding laboratory equipment.
6.2	The weight of carton was measured by the apparatus without calibration which was provided by the factory. Our inspector required the factory to calibrate the apparatus regularly.
6.3	The inspected cartons were not sealed because the factory did not repack them immediately after inspection.
6.4	Manufacturer did not allow us to draw any production sample for reference.


6.5	Our inspector required the factory to repair all defective samples before shipment.
6.6	Draft report has been provided to the factory and defective points have been explained as well.
6.7	Regarding those check items that no data or specification had been provided by the client, only measured data, which were generally judged as conformed, were listed in the data sheet for your reference.

7. Photos

	
Warehouse	Warehouse
	
Warehouse	Warehouse
	
Warehouse	Fitting test

 <p>5.1.2e</p>	
 <p>5.1.2e</p>	
	

	
<p>Fabric weight check</p>	<p>Fabric weight check</p>
	
<p>Fabric weight check</p>	<p>Bacterial filtration efficiency(BFE) test and Differential pressure test</p>
	
<p>Bacterial filtration efficiency(BFE) test and Differential pressure test</p>	<p>Bacterial filtration efficiency(BFE) test and Differential pressure test</p>

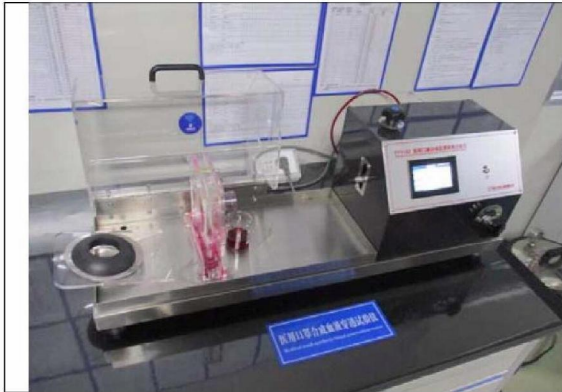


Synthetic blood penetration test




Synthetic blood penetration test



Synthetic blood penetration test

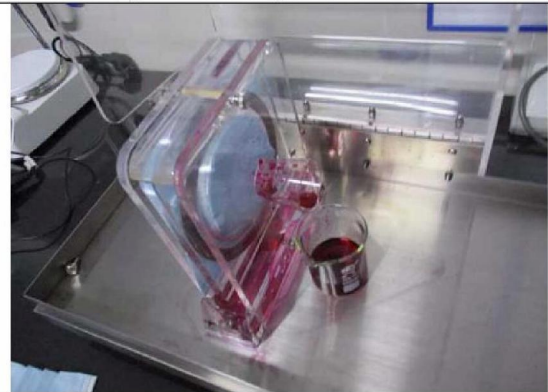

Synthetic blood penetration test

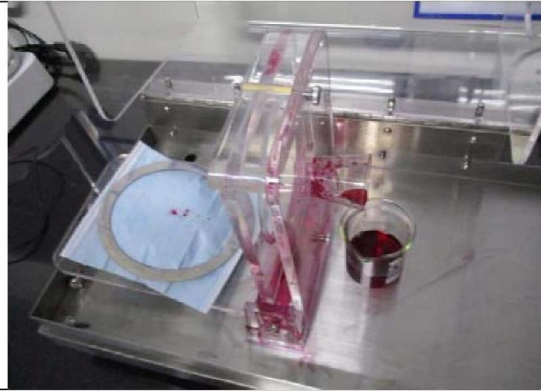


Synthetic blood penetration test




Synthetic blood penetration test




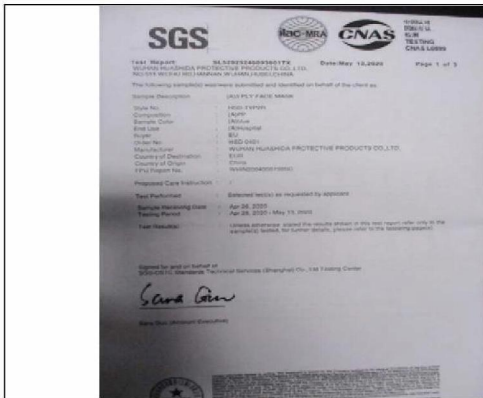

Microbial cleanliness test report

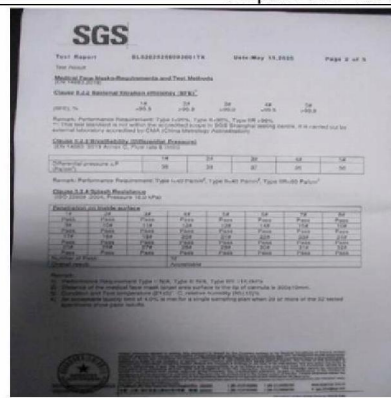

Microbial cleanliness test report

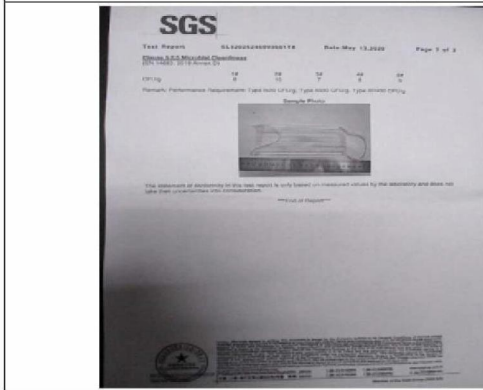

Microbial cleanliness test report



Item 4.05.053.00



Item 4.05.053.00



Item 4.05.053.00



Item 4.05.053.00

Item 4.05.053.00



Item 4.05.053.00



Item 4.05.053.00



Item 4.05.053.00




Item 4.05.053.00


<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>								
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>								
<table border="1" data-bbox="335 1142 766 1265"> <tr> <td>Manufacture date</td> <td>2020/06</td> </tr> <tr> <td>Expiry date</td> <td>2023/06</td> </tr> <tr> <td>Lot NO</td> <td>H001 20200601</td> </tr> <tr> <td>PO-NO.</td> <td>126759/LCH2020-0195-04</td> </tr> </table> <p>Manufacturer: Majestic Safety Products & Services B.V., P.O. Box B0, 3200 AB Spijkenisse, The Netherlands</p>	Manufacture date	2020/06	Expiry date	2023/06	Lot NO	H001 20200601	PO-NO.	126759/LCH2020-0195-04	<p>EN For more information, you can contact your supplier or the manufacturer. NL Voor meer informatie kunt u contact opnemen met uw leverancier of de fabrikant. FR Pour de plus amples informations vous pouvez contacter votre fournisseur ou le fabricant. DE Für weitere Informationen wenden Sie sich bitte an Ihren Lieferanten oder den Hersteller.</p>
Manufacture date	2020/06								
Expiry date	2023/06								
Lot NO	H001 20200601								
PO-NO.	126759/LCH2020-0195-04								
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>								



Item 4.05.053.00

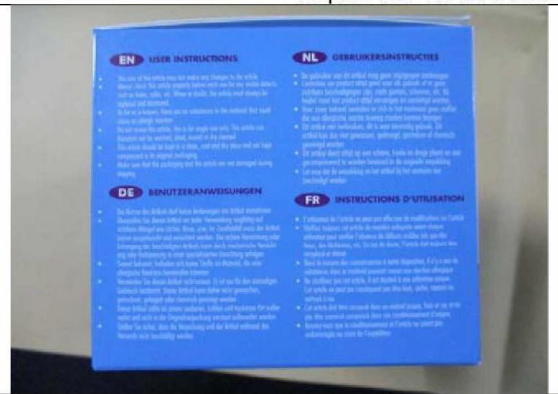


Item 4.05.053.00

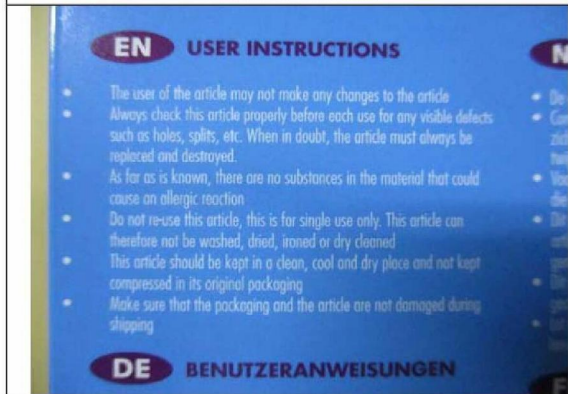


Item 4.05.053.00



Item 4.05.053.00

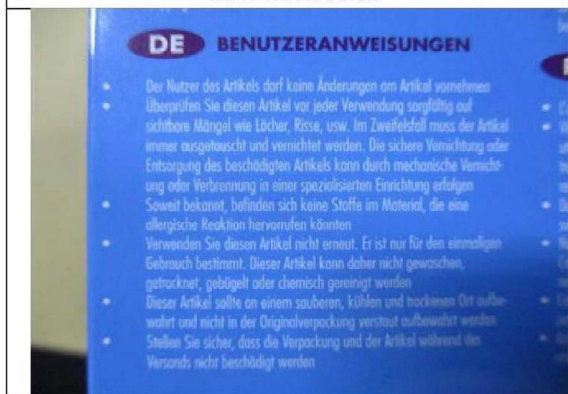

Item 4.05.053.00

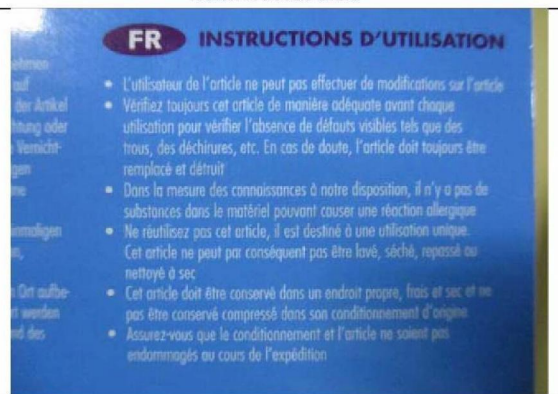

Item 4.05.053.00



Item 4.05.053.00

Item 4.05.053.00



Item 4.05.053.00



Item 4.05.053.00



Item 4.05.053.00



Item 4.05.053.00

<p>MSAFE MODEL 5053</p> <ul style="list-style-type: none">• Integrated nose bridge• Non-woven• Without glass fibres• Hypoallergenic• Very low resistance to breathing• High filtration capacity <p>Art.No. 4.05.053.00</p> <p>8 718249 066689</p>	<p>MSAFE MODEL 5053</p> <ul style="list-style-type: none">• Integrated nose bridge• Non-woven• Without glass fibres• Hypoallergenic• Very low resistance to breathing• High filtration capacity <p>Art.No. 4.05.053.00</p> <p>8 718249 066689</p>
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>
<p>Art.No. 4.05.053.00</p> <p>8 718249 066689</p>	<p>Item 4.05.053.00</p>
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>

<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>

	
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>
	
<p>Item 4.05.053.00</p>	<p>Item 4.05.053.00</p>
	<p>NA</p>
<p>Item 4.05.053.00</p>	

<p>Inspector:</p>	<p>5.1.2e</p>	<p>Auditor:</p>	<p>5.1.2e</p>
<p>End of Report</p>			