


<b>HQTS</b>	<b>INSPECTION REPORT</b> Pre-Shipment Inspection	Report No: H-20065688
-------------	---	-----------------------

INSPECTION INFORMATION			
<b>Client</b>	MEDIQ		
<b>Supplier</b>	Hangzhou Qianyuan Technology Co., Ltd		
<b>Factory</b>	Evereast Medical Products Group CO.,Ltd		
<b>Inspection Instruction</b>	GWI-H-004-HQT		
<b>P.O. No.</b>	LCH-2020-0207-05	<b>P.O. Quantity</b>	5.1.1c
<b>Item No.</b>	-	<b>Sampling Plan</b>	ISO2859 : Single, Normal
<b>Inspection Date</b>	30/06/2020 to 01/07/2020	<b>Inspection Level</b>	G-II
<b>Inspection Sequence</b>	Original inspection	<b>Sample Size</b>	800 Pieces
<b>Inspection Location</b>	Changyuan, China	<b>AQL</b>	Critical:0, Major:1.5, Minor:4.0
<b>Product Description</b>	Disposable Medical Surgery Mask		

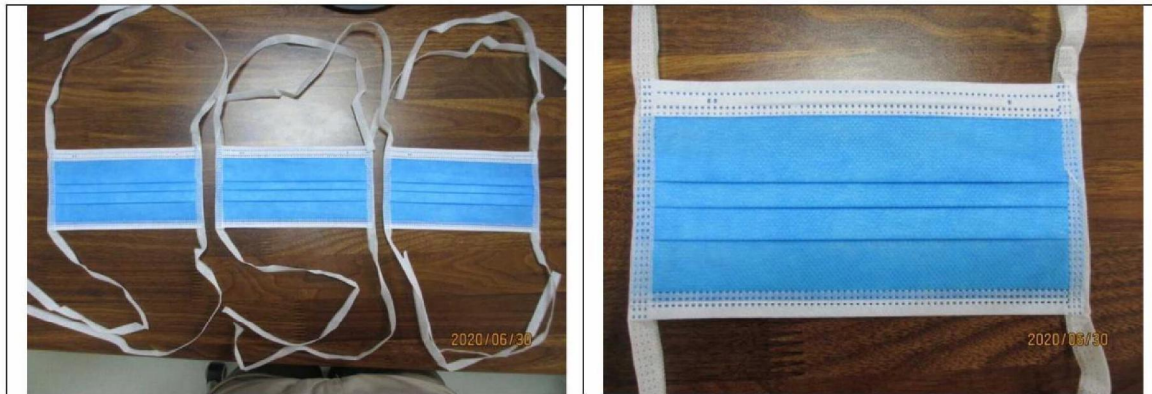

INSPECTION RESULT SUMMARY			
Category	Conform	Not Conform	Pending
1.Quantity	X		
2.Workmanship	X		
3.Style, Color & Documentation	X		
4.Field Test & Measurement			X
5.Shipping Mark & Packaging			X
<b>OVERALL CONCLUSION</b>	<b>Pending</b> for client's evaluation		

<b>Approved by HQTS Supervisor:</b>	5.1.2e
-------------------------------------	--------


*This report reflects the facts as recorded by HQTS at the time and place of inspection. It does not relieve the manufacturers from their contractual obligations nor prejudice client's right for compensation for any apparent and/or hidden defects not detected during our random inspection or occurring thereafter.*

**SPECIAL ATTENTION POINT (SAP)**

<b>1</b>	<b>Quantity</b>	NIL
<b>2</b>	<b>Workmanship</b>	NIL
<b>3</b>	<b>Style, Color &amp; Documentation</b>	NIL
<b>4</b>	<b>Field Test &amp; Measurement</b>	
4.1	Some data measurement was found out of client's specification, please refer to the data sheet marked with red.	
<b>5</b>	<b>Shipping Mark &amp; Packaging</b>	
5.1	One side of the main mark information was covered by the white label sticker, please refer to the following photos.	

	
SAP5.1-Actual found main mark information covered by white label	SAP5.1- Actual found complete main mark information

**INSPECTION DETAILS**

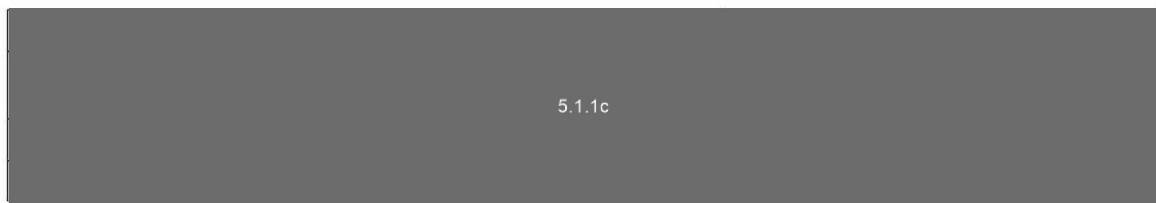

2. Workmanship	<b>Conform</b>		
Defects Description	Critical	Major	Minor
Different length of ear string(1-2.2cm)			5
Poor melting of ear string			3
Total Found:	0	0	8
Allowed:	0	21	21
Sample Size:	800	800	315

Defect Photos	
<p>Defect-Different length of ear string(1cm)-Mi</p>	<p>Defect-Different length of ear string(1.5cm)-Mi</p>
<p>Defect-Different length of ear string(2.2cm)-Mi</p>	<p>Defect-Poor melting of ear string-Mi</p>
	<p>NA</p>
<p>Defect-Poor melting of ear string-Mi</p>	

Report No: H-20065688

3. Style, Color & Documentation		Conform		
DESCRIPTION		Conform	Not Conform	Pending
The inspected samples were consistent with		-	-	-
Product Specification		X		
Digital photo		X		
Instruction manual and Document Verification		-	-	-
Client Comment		-	-	-

4. Field Test & Measurement				Pending		
Item	Specification	Sample Size	Conform	Not Conform	Pending	
4.1	Carton gross weight check	-	1 carton		X	
4.2	Carton dimension check	-	1 carton		X	
4.3	Carton drop check	-	1 carton	X		
4.4	Tension test on ear string		5 samples	X		
4.5	Product dimension check		5 samples		X	
4.6	Product weight check		5 samples	X		
4.7	Smell check		5 samples	X		
4.8	Wear check		5 samples	X		
4.9	Water leakage test		5 samples	X		
4.10	Rub test		5 samples	X		
4.11	Quantity per box check		2 Cartons	X		
4.12	Bacterial filtration efficiency test	For reference	2 samples			
4.13	NaCl particulate matter filtration efficiency test (PFE)	For reference	2 samples	X		
4.14	Synthetic blood penetration test		5 samples	X		
4.15	Differential pressure test		2 samples	X		
4.16	Microbial cleanliness test	For reference	2 samples			
4.17	Fire proof to melt-blown fabric test		5 samples	X		

Checking Item	Specification	Result
Carton gross weight check (kg) ( $\pm 10\%$ )	12	13.54
Carton dimension check (cm) ( $\pm 5\%$ )	63x45x37	63.5x46.2x39.5
Unit package weight (g)	-	99.34, 99.38, 98.25, 99.32, 98.67
Product dimension check(cm)	17.5x17	17.5x16.0/17.6x16.5 17.6x16.5/17.5x16.0 17.6x16.2
Product weight check(g)	-	4.48, 4.49, 4.49, 4.49, 4.49
Fabric weight (gsm)	-	The first layer of cloth: 39 gsm The middle layer of cloth: 28 gsm The third layer of cloth: 27 gsm

Report No: H-20065688

Differential pressure test (Pa/(cm <sup>2</sup> ))	-	54.7, 49.8
Nacl particulate matter filtration efficiency test(PFE)(0.3um)	-	96.12%, 96.24%
Test for strength of belt	>10(N)	First angle: 32N,29N Secondangle:26N,28N Third Angle:27N, 29N Fourth Angle:30N, 28N
Microbial cleanliness test	≤30cfu/g	
Bacterial filtration efficiency test		

<b>5. Shipping Mark &amp; Packaging</b>	<b>Pending</b>
---	----------------

5.1 Shipping mark photos



5.2 Packaging photos



Packaging Description:

- 20 pieces of products were packed into an inner bag.
- 125 bags were packed into a 5-ply corrugated export carton.
- Export carton was sealed with gummed tape and secured with nylon strap.

<b>6. Reference Note</b>	
6.1	For bacterial filtration efficiency test, The test process was not supervised because it needed 48 hours of bacterial and microbial culture, While Nacl particulate matter filtration efficiency test(PFE) was conducted on site for your reference. And our inspector sampling to factory, the test process was not supervised by our inspector and the result would be provided later by the factory.
6.2	For microbial cleanliness test, The test process was not supervised because it needed 48 hours of bacterial and microbial culture, and our inspector sampling to factory, the test process was not supervised by our inspector and the result would be provided later by the factory.
6.3	Manufacturer did not allow us to draw any production sample for reference.
6.4	Our inspector required the factory to repair all defective samples before shipment.
6.5	For bacterial filtration efficiency test & microbial cleanliness test, the factory claimed those test will take about 3 days, and the test result will be available on 3 <sup>rd</sup> .Jul.2020.

7. Photos



Warehouse



Drawn cartons




Carton gross weight check



Carton size check



Product size check

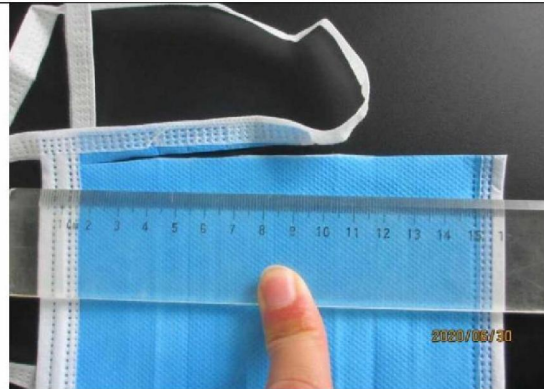

Product size check



Differential pressure test & NaCl particulate matter filtration efficiency test

Differential pressure test & NaCl particulate matter filtration efficiency test

Qualification certificate


Synthetic blood penetration test

Synthetic blood penetration test

Synthetic blood penetration test



Qualification certificate

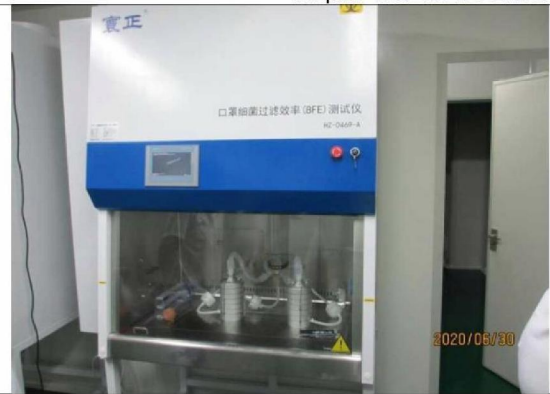

Bacterial filtration efficiency test



Bacterial filtration efficiency test



Bacterial filtration efficiency test



Qualification certificate



Microbial cleanliness test


 <p>注册地址：济南市天桥区蓝翔中路5号天梯大厦 Calibration 校准证 pass 器具编号S/N: BSC11B15070460 校准日期Cal.Date: 2019年9月24日 建议复校日期DUE: 2020年9月23日 校准员Cal.by: 10 校准单位：深圳新广行检测技术有限公司 电话：0755-26450356 2020/06/30</p>	 <p>METTLER TOLEDO 44897.0g 2020/06/30</p>
<p>Qualification certificate</p>	<p>Product weight check</p>
 <p>METTLER TOLEDO 0.2853g 2020/06/30</p>	 <p>METTLER TOLEDO 0.3997g 2020/06/30</p>
<p>Fabric weight check</p>	<p>Fabric weight check</p>
 <p>METTLER TOLEDO 0.2719g 2020/06/30</p>	 <p>Calibration 校准证 pass 编号S/N: M11785758767 校准日期Cal.Date: 2019年9月26日 建议复校日期DUE: 2020年9月23日 校准单位：深圳新广行检测技术有限公司 电话：0755-26450356 2020/06/30</p>
<p>Fabric weight check</p>	<p>Calibration report</p>



 <p>2020/07/01</p>	 <p>2020/06/30</p>
<p>Rub test</p>	<p>Factory</p>
 <p>2020/06/30</p>	 <p>2020/06/30</p>
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>
 <p>2020/06/30</p>	 <p>2020/06/30</p>
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>

 <p>LCH-2020-0207-05          Order number:2020024-5          Mediq Nederland B.V.          Acting as part of: the Landelijk Consortium Hulpmiddelen (LCID)          Rijzatho 10          3454 PE Utrecht          The Netherlands          Disposable Sterile Surgical Mask          CARTON NO. 1501-2300          MADE IN CHINA</p> <p style="text-align: right; color: orange;">2020/06/30</p>	 <p>LOT: 200605          MFG: 2020-06-05          EXP: 2023-06-04          G.W.: 12kg          N.W.: 11kg          MEAS: 63×45×37cm</p> <p>↑↑    ☉    ☀    ☔</p> <p style="text-align: right; color: orange;">2020/06/30</p>
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>
 <p>EN14683-2019 TYPE IIR          BFE≥98%</p> <p><b>Disposable Sterile Surgical Mask</b>          Made in China</p> <p>2500PCS</p> <p>Shanghai International Holding Corp. GmbH (Europe)          Add: Eiffelstrasse 80, 20537 Hamburg, Germany</p> <p style="text-align: right; color: orange;">2020/06/30</p>	 <p style="text-align: right; color: orange;">2020/06/30</p>
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>
 <p style="text-align: right; color: orange;">2020/06/30</p>	 <p style="text-align: right; color: orange;">2020/06/30</p>
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>

	
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>
	
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>
	
<p>Item Disposable medical surgery mask</p>	<p>Item Disposable medical surgery mask</p>

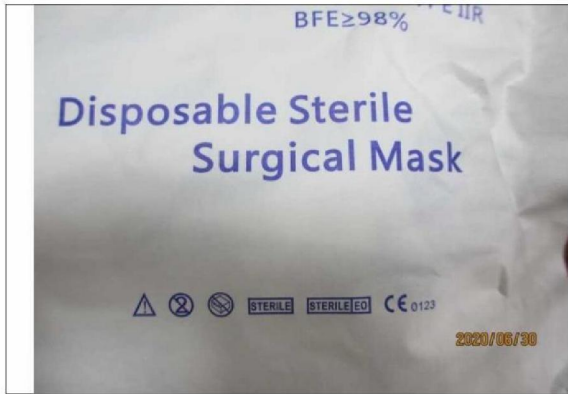

Item Disposable medical surgery mask



Item Disposable medical surgery mask



Item Disposable medical surgery mask




Item Disposable medical surgery mask




Item Disposable medical surgery mask



Item Disposable medical surgery mask

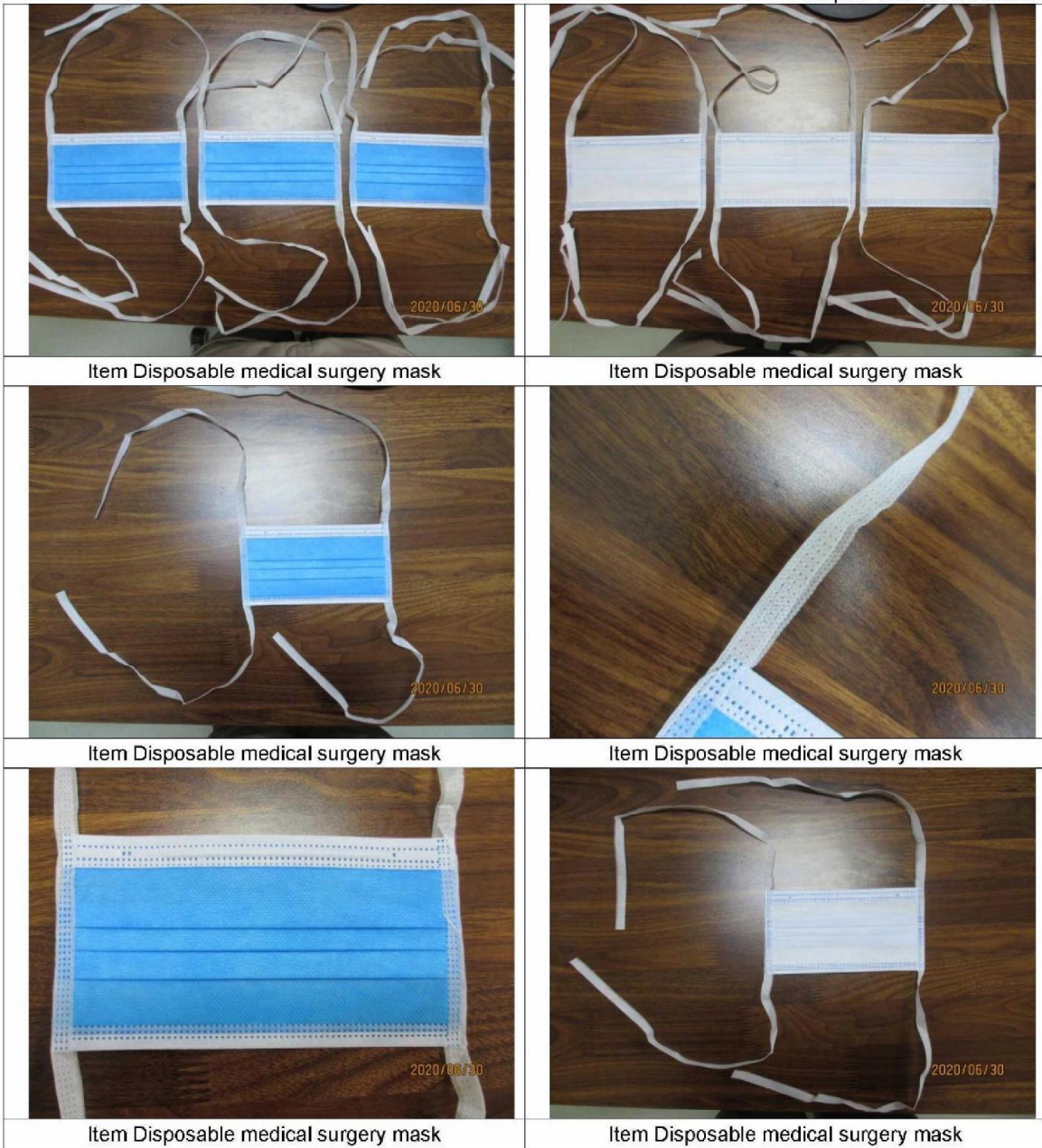

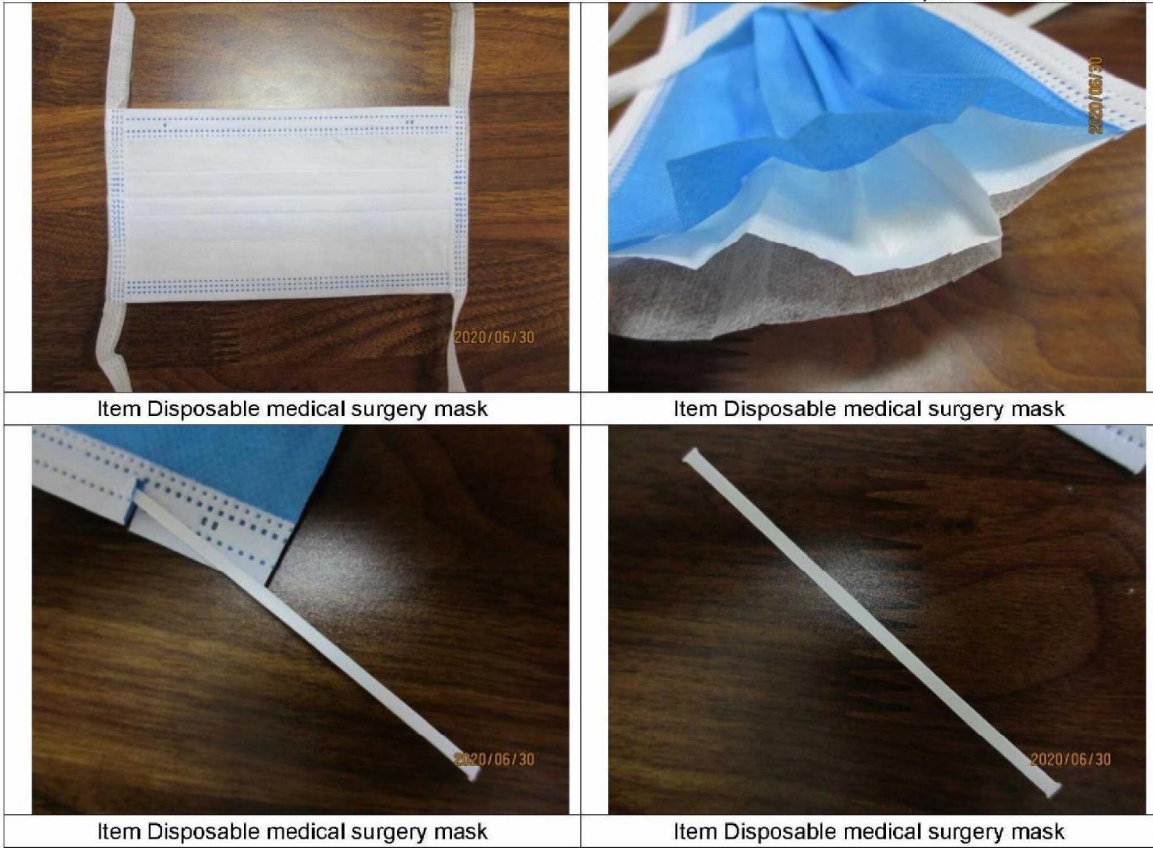

<b>Inspector:</b>	5.1.2e	<b>Auditor:</b>	5.1.2e
<b>End of Report</b>			