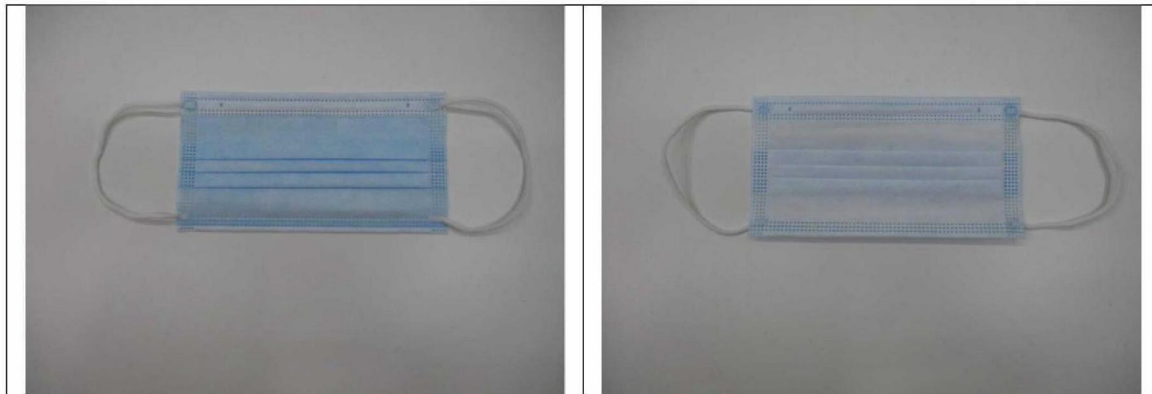


HQTS	INSPECTION REPORT Pre-Shipment Inspection	Report No: H-20077290
-------------	-----------------------------------------------------	-----------------------

INSPECTION INFORMATION			
Client	MEDIQ		
Supplier	Bingo Times Guangzhou Supply Chain Co., Ltd		
Factory	Bingo Times Guangzhou Supply Chain Co., Ltd		
Inspection Instruction	GWI-H-004-HQT&SWI-2099-HQT-medi		
P.O. No.	LCH2020-0222-05	P.O. Quantity	5.1.1c
Item No.	KZ-05	Sampling Plan	ISO2859 : Single, Normal
Inspection Date	13/07/2020	Inspection Level	G-II
Inspection Sequence	Original inspection	Sample Size	800 Pieces
Inspection Location	Guangzhou, China	AQL	Critical:0, Major:1.5, Minor:4.0
Product Description	SURGICAL MASK TYPE IIR		




INSPECTION RESULT SUMMARY			
Category	Conform	Not Conform	Pending
1.Quantity	X		
2.Workmanship			X
3.Style, Color & Documentation			X
4.Field Test & Measurement	X		
5.Shipping Mark & Packaging	X		
OVERALL CONCLUSION	PENDING for client's evaluation		

Approved by HQTS Supervisor:	5.1.2e
-------------------------------------	--------


This report reflects the facts as recorded by HQTS at the time and place of inspection. It does not relieve the manufacturers from their contractual obligations nor prejudice client's right for compensation for any apparent and/or hidden defects not detected during our random inspection or occurring thereafter.

SPECIAL ATTENTION POINT (SAP)

1	Quantity	NIL
2	Workmanship	
2.1	Actual found two production processes were found in inspected samples, lot A was about 85%, lot B was about 15%, for details, please refer to photo attached.	

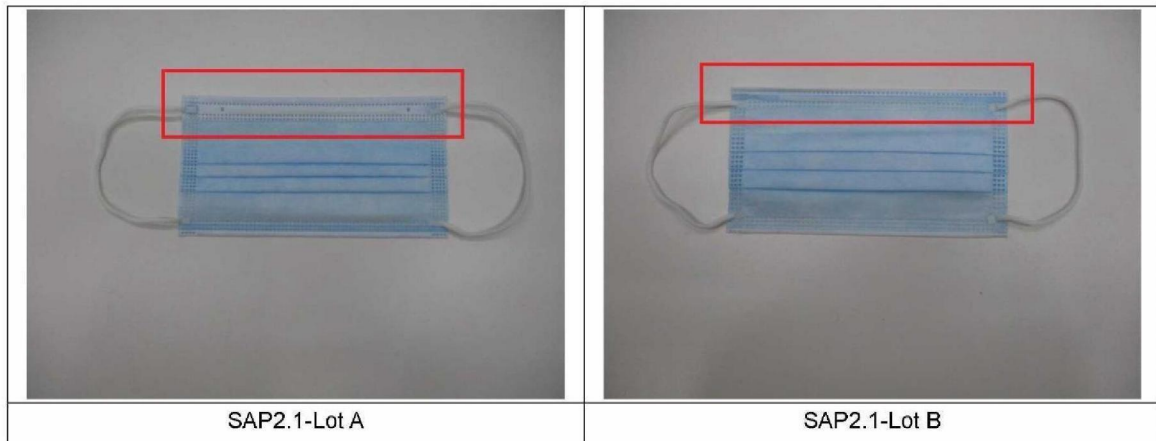

3	Style, Color & Documentation	
3.1	Actual found the opening direction of color box was difference form spec, for details, please refer to photo attached.	
3.2	Actual found the content printed on the color box was different form spec, for details, please refer to photo attached.	

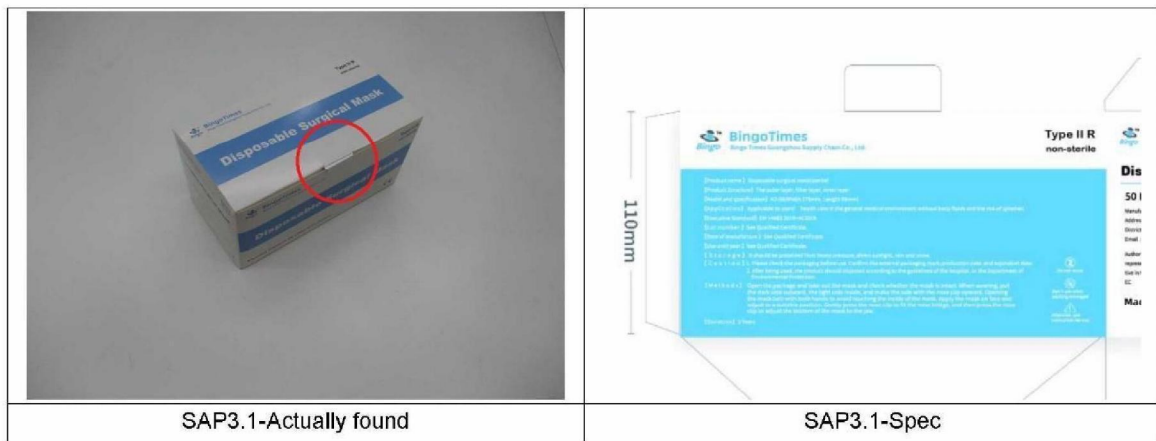


<p> [Product name] Disposable surgical mask(non-sterile) [Product Structure] The outer layer, filter layer, inner layer [Model and specification] KZ-05(Width 175mm, Length 95mm) [Application] Applicable to users' health care in the general medical env [Executive Standard] EN 14683 2019+AC2019 [Lot number] See Package KZ0520200705 [Date of manufacture] See Package 2020-07-05 [Use until year] See Package 2022-07-04 [Storage] It should be protected from heavy pressure, direct sunlight, r [Caution] 1. Please check the packaging before use. Confirm the external 2. After being used, the product should disposed according to t Environmental Protection. [Methods] Open the package and take out the mask and check wheti the dark side outward, the light side inside, and make the the mask held with both hands to avoid touching the inside </p>	<p> [Product name] Disposable surgical mask(sterile) [Product Structure] The outer layer, filter layer, inner layer [Model and specification] KZ-05(Width 175mm, Length 95mm) [Application] Applicable to users' health care in the general media [Executive Standard] EN 14683 2019+AC2019 [Lot number] See Qualified Certificate. [Date of manufacture] See Qualified Certificate. [Use until year] See Qualified Certificate. [Storage] It should be protected from heavy pressure, direct su [Caution] 1. Please check the packaging before use. Confirm the </p>
<p>SAP3.2-Actually found</p>	<p>SAP3.2-Spec</p>

<p>4</p>	<p>Field Test & Measurement</p>	<p>NIL</p>
----------	-------------------------------------	------------

<p>5</p>	<p>Shipping Mark & Packaging</p>	<p>NIL</p>
----------	--------------------------------------	------------

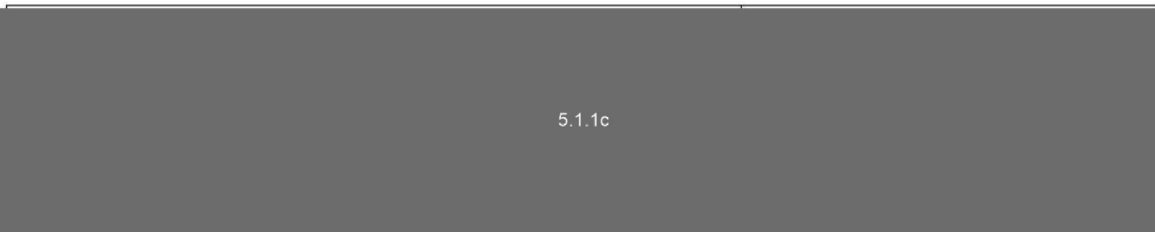INSPECTION DETAILS



5.1.1c


2. Workmanship	Pending		
Defects Description	Critical	Major	Minor
Asymmetrical of ear loop welding dot(3mm)		5	
Pull apart of belt(use pull with normal strength)		2	
Dirty mark on belt		2	
Short-long belt(1.5cm)		1	
Folding mark on inner(4-5cm)			2
Impurity mark on surface(1mm)			2
Slightly wrinkle on inner(1.2cm)			1
Total Found:	0	10	5
Allowed:	0	21	21
Sample Size:	0	800	315

<p>SPECIAL ATTENTION POINT:</p>	
<p>In addition to defects recorded in the workmanship table, the below massive defects were also found:</p>	
<p>2.1</p>	<p>Actual found two production processes were found in inspected samples, lot A was about 85%, lot B was about 15%</p>

Defect Photos	
	
<p>Defect-Asymmetrical of ear loop welding dot(3mm)-ma (2)</p>	<p>Defect-Asymmetrical of ear loop welding dot(3mm)-ma (3)</p>
	
<p>Defect-Asymmetrical of ear loop welding dot(3mm)-ma (1)</p>	<p>Defect-Pull apart of belt(use pull with normal strength)-ma</p>
	
<p>Defect-Pull apart of belt(use pull with normal strength)-ma</p>	<p>Defect-dirty mark on belt-ma (1)</p>

<p>Defect-dirty mark on belt-ma (2)</p>	<p>Defect-short-long belt(1.5cm)-ma</p>
<p>Defect-folding mark on inner(4-5cm)-mi (1)</p>	<p>Defect-folding mark on inner(4-5cm)-mi (2)</p>
<p>Defect-impurity mark on surface(1mm)-mi (1)</p>	<p>Defect-impurity mark on surface(1mm)-mi (2)</p>

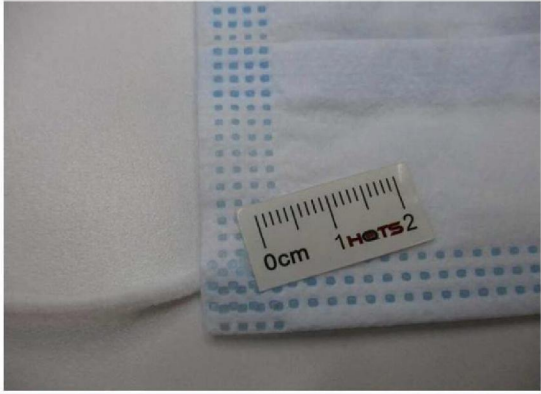Report No: H-20077290



	NA
Defect-slightly wrinkle on inner(1.2cm)-mi	

3. Style, Color & Documentation	Pending		
DESCRIPTION	Conform	Not Conform	Pending
The inspected samples were consistent with	-	-	-
Purchase Order	X		
Digital photo			X
Client Comment	-	-	-

4. Field Test & Measurement			Conform		
Item	Specification	Sample Size	Conform	Not Conform	Pending
4.1	Carton dimension check	1 carton	X		
4.2	Carton gross weight check	1 carton	X		
4.3	Carton drop check	1 carton	X		
4.4	Quantity per carton check	All selected cartons	X		
4.5	Product dimension check	5 samples	X		
4.6	Product net weight check	5 samples	X		
4.7	Smell check	5 samples	X		
4.8	Water leakage test	5 samples	X		
4.9	Wear check	5 samples	X		
4.10	Rub test	5 samples	X		
4.11	Tension test on ear string	1kg	X		
4.12	Bacterial filtration efficiency test	5 samples	-	-	-
4.13	Synthetic blood penetration test	5 samples	-	-	-
4.14	Differential pressure test	5 samples	X		
4.15	Microbial cleanliness test	5 samples	-	-	-
4.16	Fire proof to melt-blown fabric test	5 samples	X		
4.17	Nacl particulate matter filtration efficiency test (PFE)	0.3um	X		

Report No: H-20077290

Checking Item	Specification	Result
Carton dimension (mm) ($\pm 5\%$)	545x460x415	543x421x465
Carton gross weight (kg)/ ($\pm 10\%$)	10.08	9.672
Color box dimension(mm)	200x105x110	200x105x111, 200x105x111, 200x105x111, 200x105x111, 200x105x111
Color box weight (g)	-	214.6, 217.8, 216.3, 215.9, 218.2
Product dimension(mm)	175x95	175x95, 175x95, 175x95, 175x95, 175x95
Product net weight (g)	-	3.5, 3.6, 3.5, 3.5, 3.5
Fabric weight (gsm)	The first layer of cloth	25.81, 25.32, 24.85, 24.91, 25.14
	The melt-blown fabric	25.27, 25.12, 25.32, 24.89, 25.18
	The third layer of cloth	24.68, 25.19, 24.93, 24.55, 25.24
Differential pressure test (Pa/(cm ²))	-	34.86, 35.12, 34.97, 35.23, 34.89
Particulate matter filtration efficiency test(PFE)(0.3um)	-	98.05%, 97.53%, 99.01%, 98.06%, 98.39%

5. Shipping Mark & Packaging	Conform
5.1 Shipping mark photos	
	
5.2 Packaging photos	

Report No: H-20077290

**Packaging Description:**

- 25 pieces of products were packed into a polybag.
- 2 polybags were packed into a gift box.
- 40 gift boxes were packed into a 5-ply corrugated export carton.
- Export carton was sealed with gummed tape.

6. Reference Note

6.1	The weight was measured by the apparatus without calibration which was provided by the factory. Our inspector required the factory to calibrate the apparatus regularly.
6.2	The factory could not provide relevant apparatus so the Bacterial filtration efficiency test, Synthetic blood penetration test and Microbial cleanliness test were not conducted during the inspection.
6.3	Carton number was not assigned by the factory yet.
6.4	Actual found the export carton size printing on shipping was LXHXW instead of LXWXH.
6.5	100 pieces of normal samples were drawn and kept in the factory due to some tests were not done during inspection.
6.6	The factory claimed that all inspected sample(800pieces) would be destroyed.

7. Photos



Warehouse



Warehouse



Warehouse



Warehouse



Warehouse



Item KZ-05

<p>Item KZ-05</p>	<p>Item KZ-05</p>
<p>Item KZ-05</p>	<p>Item KZ-05</p>
<p>Item KZ-05</p>	<p>Item KZ-05</p>

<p>Instruction for User [Product name] Disposable surgical mask(non-sterile) [Application] Applicable to users' health care in the general medical environment without body fluids and the risk of splashes. [Methods] Open the package and take out the mask and check whether the mask is intact. When wearing, put the dark side outward. Adjust the nose bridge, and make the side with the nose clip upward. Open the mask belt with both hands to avoid touching the inside of the mask. Apply the mask on face and adjust it to a suitable position. Gently press the nose clip to fit the nose bridge, and then press the nose clip to adjust the bottom of the mask to the jaw. [Duration] 2 years [Storage] It should be protected from heavy pressure, direct sunlight, rain and snow. [Note] 1. Please choose sterile or non-sterile products according to requirements. 2. This product is disposable. Second use is strictly prohibited. It should be disposed of in accordance with the requirements of hospitals and environmental protection departments when used. 3. Please check whether the packaging is intact before use. Do not use the product if the package or mask is damaged. Confirm the external packaging mark, production date and validity period. Use within the validity period. 4. Be caution using the product if you are allergic to non-woven fabrics. 5. During use, the product should be replaced in time in the following situations: a significant increase in respiratory resistance (when breathing is difficult), the broken mask, and the contaminative mask (when contaminated with foreign matters such as blood stains and droplets). Manufacturers: Bingo Times Guangzhou Supply Chain Co.,Ltd Address: 101 No. 13 Middle Daguan Road Tianhe District Guangzhou Guangdong China Email: 5.12e@bingostyle.cn Tel: 020-35632599 CompanySUNGO Europe B.V. AddressOlympisch Stadion 24, 1076DE Amsterdam, Netherlands Tel: +31 (0) 20211106 E-mail: 5.12e@sungogroup.com Don't use when packing damaged Date: Apr.22,2020 Version:A/0</p>	<p>Instruction for User [Product name] Disposable surgical mask(non-sterile) [Application] Applicable to users' health care in the general medical environment without body fluids and the risk of splashes. [Methods] Open the package and take out the mask and check whether the mask is intact. When wearing, put the dark side outward, the light side inside, and make the side with the nose clip upward. Opening the mask belt with both hands to avoid touching the inside of the mask. Apply the mask on face and adjust to a suitable position. Gently press the nose clip to fit the nose bridge, and then press the nose clip to adjust the bottom of the mask to the jaw. [Duration] 2 years [Storage] It should be protected from heavy pressure, direct sunlight, rain and snow. [Note] 1. Please choose sterile or non-sterile products according to requirements. 2. This product is disposable. Second use is strictly prohibited. It should be disposed of in accordance with the requirements of hospitals and environmental protection departments when used. 3. Please check whether the packaging is intact before use. Do not use the product if the package or mask is damaged. Confirm the external packaging mark, production date and validity period. Use within the validity period. 4. Be caution using the product if you are allergic to non-woven fabrics.</p>
<p>Item KZ-05</p> <p>[Duration] 2 years [Storage] It should be protected from heavy pressure, direct sunlight, rain and snow. [Note] 1. Please choose sterile or non-sterile products according to requirements. 2. This product is disposable. Second use is strictly prohibited. It should be disposed of in accordance with the requirements of hospitals and environmental protection departments when used. 3. Please check whether the packaging is intact before use. Do not use the product if the package or mask is damaged. Confirm the external packaging mark, production date and validity period. Use within the validity period. 4. Be caution using the product if you are allergic to non-woven fabrics. 5. During use, the product should be replaced in time in the following situations: a significant increase in respiratory resistance (when breathing is difficult), the broken mask, and the contaminative mask (when contaminated with foreign matters such as blood stains and droplets). Manufacturers: Bingo Times Guangzhou Supply Chain Co.,Ltd Address: 101 No. 13 Middle Daguan Road Tianhe District Guangzhou Guangdong China Email: 5.12e@bingostyle.cn Tel: 020-35632599 CompanySUNGO Europe B.V. AddressOlympisch Stadion 24, 1076DE Amsterdam, Netherlands Tel: +31 (0) 20211106 E-mail: 5.12e@sungogroup.com Don't use when packing damaged Date: Apr.22,2020 Version:A/0</p>	<p>Item KZ-05</p> <p>5. During use, the product should be replaced in time in the following situations: a significant increase in respiratory resistance (when breathing is difficult), the broken mask, and the contaminative mask (when contaminated with foreign matters such as blood stains and droplets). Manufacturers: Bingo Times Guangzhou Supply Chain Co.,Ltd Address: 101 No. 13 Middle Daguan Road Tianhe District Guangzhou Guangdong China Email: 5.12e@bingostyle.cn Tel: 020-35632599 CompanySUNGO Europe B.V. AddressOlympisch Stadion 24, 1076DE Amsterdam, Netherlands Tel: +31 (0) 20211106 E-mail: 5.12e@sungogroup.com Don't use when packing damaged Date: Apr.22,2020 Version:A/0</p>
<p>Item KZ-05</p> <p>Disposable Surgical Mask</p>	<p>Item KZ-05</p> <p>BingoTimes Disposable Surgical Mask Type II R non-sterile</p>

<p>Item KZ-05</p>	<p>Item KZ-05</p>
<p>Item KZ-05</p>	<p>Item KZ-05</p>
<p>Item KZ-05</p>	<p>Item KZ-05</p>

<p>[Application] Applicable to users' health care in the general medical environment.</p> <p>[Executive Standard] EN 14683 2019+AC2019</p> <p>[Lot number] See Package KZ0520200705</p> <p>[Date of manufacture] See Package 2020-07-05</p> <p>[Use until year] See Package 2022-07-04</p> <p>[Storage] It should be protected from heavy pressure, direct sunlight.</p> <p>[Caution] 1. Please check the packaging before use. Confirm the external packaging mark, production date, and expiration date. 2. After being used, the product should be disposed according to Environmental Protection.</p> <p>[Methods] Open the package and take out the mask and check whether the mask is intact. When wearing, put the dark side outward, the light side inside, and make the mask belt with both hands to avoid touching the inside of the mask. Gently press the nose clip to adjust the bottom of the mask to the jaw.</p> <p>[Duration] 2 Years</p>	<p>(mm)</p> <p>general medical environment without body fluids and the risk of splashes.</p> <p>KZ0520200705</p> <p>2020-07-05</p> <p>2022-07-04</p> <p>direct sunlight, rain and snow.</p> <p>Confirm the external packaging mark, production date, and expiration date. Dispose according to the guidelines of the hospital, or the Department of Environmental Protection.</p> <p>mask and check whether the mask is intact. When wearing, put the dark side outward, the light side inside, and make the mask belt with both hands to avoid touching the inside of the mask. Gently press the nose clip to fit the nose bridge, and then press the nose clip to adjust the bottom of the mask to the jaw.</p> <p>Do not reuse</p> <p>Don't use when packing damaged</p> <p>Attention, use instructions for use</p>
-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------

Item KZ-05

Item KZ-05



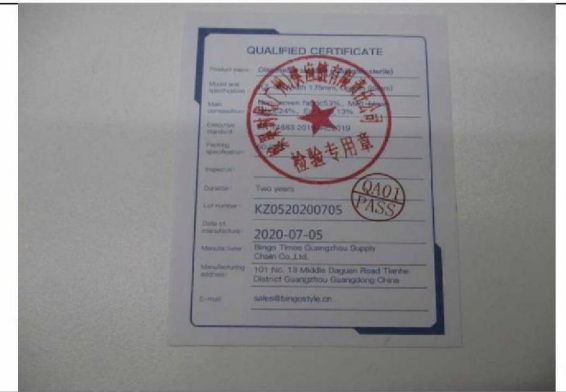
Item KZ-05




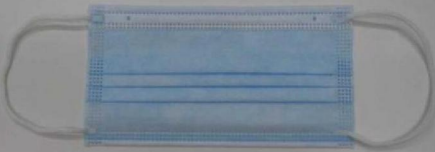




Item KZ-05



Item KZ-05



Item KZ-05


	
<p>Item KZ-05</p>	<p>Item KZ-05</p>
	
<p>Item KZ-05</p>	<p>Tension test on ear string</p>
	
<p>Water leakage test</p>	<p>Fire proof to melt-blown fabric test</p>

	
<p>Fire proof to melt-blown fabric test</p>	<p>Fabric weight check</p>
	
<p>Fabric weight check</p>	<p>Fabric weight check</p>
	
<p>Differential pressure test</p>	<p>Differential pressure test</p>




Differential pressure test

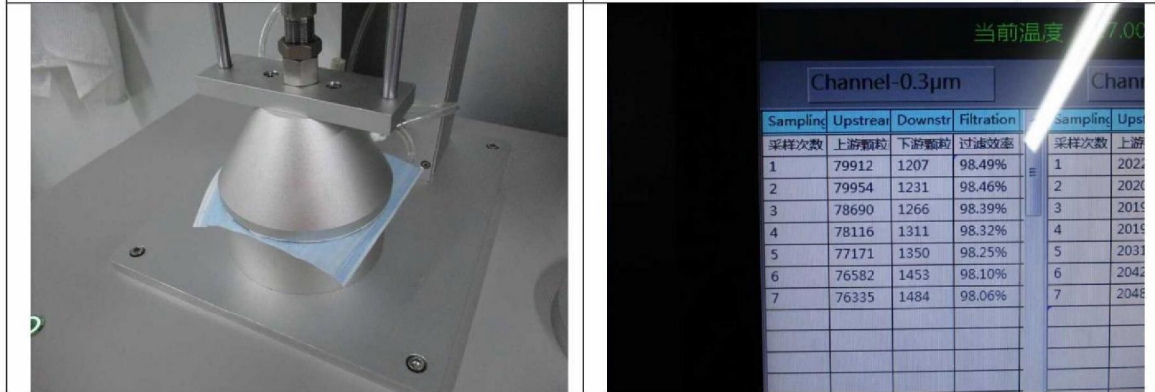NaCl particulate matter filtration efficiency test (PFE)



NaCl particulate matter filtration efficiency test (PFE)

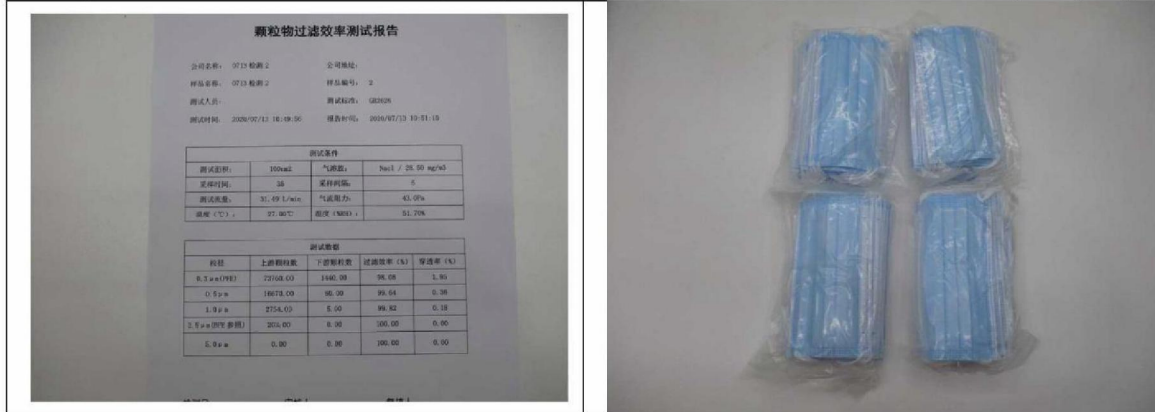NaCl particulate matter filtration efficiency test (PFE)



NaCl particulate matter filtration efficiency test (PFE)

Sealed samples

Report No: H-20077290



Inspector:	5.1.2e	Auditor:	5.1.2e
End of Report			