



In Vitro Diagnostic Reagent

Clinical Evaluation Report

Product Name: Rapid SARS-CoV-2 Antigen Test Card

Manufacturer: MP Biomedicals Germany GmbH

Report Date: February 14, 2021



Table of Contents

ABSTRACT	I
1 OVERVIEW.....	1
1.1 INTRODUCTION.....	1
1.1.1 <i>Background for Clinical Trial</i>	1
1.1.2 <i>Product Testing Principle and Methods</i>	1
1.1.3 <i>Intended Use</i>	1
1.2 PURPOSE	2
1.3 CLINICAL TRIAL MANAGEMENT.....	2
1.4 EXPERIMENT DESIGN.....	2
1.4.1 <i>Description of Overall Experiment Design and Plan</i>	2
1.4.2 <i>Sample Selection and Exclusion Criteria</i>	3
1.4.3 <i>Testing Reagent</i>	4
1.4.4 <i>Quality Control Methods</i>	4
1.4.5 <i>Statistical Analysis Methods</i>	4
1.5 CLINICAL TRIAL RESULTS AND ANALYSIS.....	5
1.5.1 <i>Sample distribution statistics</i>	5
1.5.2 <i>Statistical analysis of datas</i>	5
2 APPENDIX.....	9



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

Abstract

Chongqing Public Health Medical Center, Ankara Dr. Nafiz Korez Sincan State Hospital, and the Catholic University of Korea Uijeongbu St. Mary's Hospital carried out the clinical trials of the Rapid SARS-CoV-2 Antigen Test Card ("Evaluated Reagent").

The clinical diagnosis ability and the in-use safety and effectiveness of the Rapid SARS-CoV-2 Antigen Test Card were evaluated based on comparison with RT-PCR method.

The clinical trials were carried out in strict accordance with the experiment plan. A total of 1027 samples were used, including 313 positive samples and 714 negative samples. The test results are summarized below:

Evaluated Reagent Results	RT-PCR Results		Total
	Positive (+)	Negative (-)	
Positive (+)	301	6	307
Negative (-)	12	708	720
Total	313	714	1027

Sensitivity (PPA) = $301/313 \times 100\% = 96.17\%$ (95% CI: 94.04% - 98.29%)

Specificity (NPA) = $708/714 \times 100\% = 99.16\%$ (95% CI: 98.49% - 99.83%)

Accuracy (OPA) = $1009/1027 \times 100\% = 98.25\%$ (95% CI: 97.44% - 99.05%)

PPV = $301/307 \times 100\% = 98.05\%$ (95% CI: 96.50% - 99.59%)

NPV = $708/720 \times 100\% = 98.33\%$ (95% CI: 97.40% - 99.27%)

Detection rate of positive samples with different onset days are summarized below:

	Number of Samples	Detection Rate
1-3 days	165	161/165, 97.58%
4-7 days	148	140/148, 94.59%

Detection rate of 1-3 days = $161/165 \times 100\% = 97.58\%$ (95% CI: 95.23% - 99.92%)

Detection rate of 4-7 days = $140/148 \times 100\% = 94.59\%$ (95% CI: 90.95% - 98.24%)



MP Biomedicals Germany GmbH
Thueringer Straße 15
37269 Eschwege GERMANY
mpbio.com

1 Overview

1.1 Introduction

1.1.1 Background for Clinical Trial

COVID-19 is caused by infection with SARS-CoV-2 (novel coronavirus) with major clinical symptoms of fever, cough, and muscle fatigue. The incubation period is 1 to 14 days, mostly 3 to 7 days. A small number of patients can also exhibit nasal congestion, tears, sore throat, myalgia, diarrhea, and other symptoms. Critically ill patients often develop dyspnea and/or hypoxemia one week after disease onset, with severe cases rapidly progressing into acute respiratory distress syndrome, septic shock, metabolic acidosis, coagulation disorder, and multiple organ failure. SARS-CoV-2 is an RNA virus belonging to the β -CoV lineage. It has a viral envelope, is round or oval in shape with approximately 60-140 nm in diameter, and is usually polymorphic.

1.1.2 Product Testing Principle and Methods

1. Testing Principle

Rapid SARS-CoV-2 Antigen Test Card is an immunochromatographic lateral flow device that employs the principle of double antibody sandwich method. Colloidal gold conjugated anti-SARS-CoV-2 antibodies are dry-immobilized on the test device. When the specimen is added, it migrates by capillary diffusion through the strip to re-hydrate the gold conjugate complexes. If present at or above the limit of detection, SARS-CoV-2 viral antigens will react with the gold conjugate complexes to form particles, which will continue to migrate along the strip until the Test Zone (T) where they are captured by the immobilized anti-SARS-CoV-2 antibodies to form a visible red line. If there are no SARS-CoV-2 viral antigens in the specimen, no red line will appear in the Test Zone (T). The gold conjugate complexes will continue to migrate alone until being captured by immobilized antibody in the Control Zone (C) to form a red line, which indicates the validity of the test.

2. Methods

Read the instructions for use carefully before testing. Bring the kit components to room temperature and carry out the test under room temperature.

Refer to the instructions for use for testing procedures.

1.1.3 Intended Use

Rapid SARS-CoV-2 Antigen Test Card is an immunochromatography based one step in vitro test. It is designed for the rapid qualitative determination of SARS-CoV-2 virus antigen in nasal swabs, nasopharyngeal swabs or oropharyngeal swabs from individuals suspected of COVID-19 by their healthcare provider within the first seven days of symptom onset.

Rapid SARS-CoV-2 Antigen Test Card cannot be used as the basis to diagnose or exclude SARS-CoV-2 infection.



MP Biomedicals Germany GmbH
Thueringer Straße 15
37269 Eschwege GERMANY
mpbio.com

1.2 Purpose

To evaluate the clinical diagnosis ability and the in-use safety and effectiveness of the Rapid SARS-CoV-2 Antigen Test Card manufactured by MP Biomedicals Germany GmbH.

1.3 Clinical Trial Management

1. Quality control of clinical test conditions and personnel: Before the clinical test is carried out, the department undertaking the clinical test shall check and confirm that the test equipment and instruments are running normally, the reagent is within the effective period and the storage conditions meet the requirements of the kit. Test personnel are trained and able to meet test requirements.
2. Quality control during the clinical trial: This clinical trial adopted a blind controlled trial design, and the samples were blinded. Neither the sample tester nor the result observer knew about the blindness of the samples. It can effectively reduce the bias of test results and ensure the quality of clinical trials. The acquisition, preservation and test operation of samples shall be carried out in strict accordance with the program, and clinical supervisors shall be assigned to supervise the implementation of clinical trials.
3. Quality control of clinical trial results: control products are set to ensure the reliable quality of experimental results.
4. The main researcher should appoint experienced personnel in his department to carry out the experiment and judge the results.
5. The operation of clinical trials shall be carried out in strict accordance with standard operating procedures.

1.4 Experiment Design


1.4.1 Description of Overall Experiment Design and Plan

1. Methods

In this clinical trial, blind data analysis was adopted and the Rapid SARS-CoV-2 Antigen Test Card was used to compare and study the detection results of the samples with the RT-PCR method. The specific test process is as follows:



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

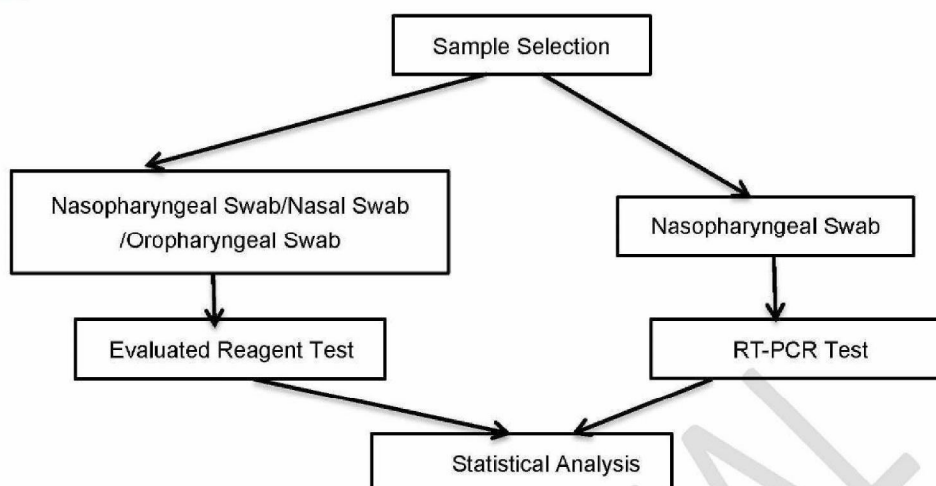

2. Selection of Clinical Trial Institutions

The clinical trials were conducted in Chongqing Public Health Medical Center, Ankara Dr. Nafiz Korez Sincan State Hospital, and the Catholic University of Korea Uijeongbu St. Mary's Hospital. The hospitals have the technical personnel and relevant experience in clinical trial projects to ensure the successful implementation of the clinical trials.

3. Number of samples

Positive samples: at least 30 cases.

Negative samples: at least 50 cases.

4. Sample type

Nasopharyngeal swab/nasal swab/oropharyngeal swab samples from the patients.

1.4.2 Sample Selection and Exclusion Criteria

1. Inclusion Criteria:

- a. COVID-19 suspected cases (confirmed cases, excluded cases) were selected.
- b. No limitation on gender and age.
- c. The patient information and clinical test/diagnosis information of the samples can be completely traced.

2. Exclusion Criteria:

- a. Insufficient sample size.
- b. Improper storage of samples.
- c. Others deemed reasonable by the principal investigator.

3. Rejection Criteria:

- a. Incomplete case history or failure to meet the inclusion criteria.



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

- b. The normal test cannot be completed due to equipment or operation factors (the sample is contaminated during operation and the sample size is insufficient).
- c. Samples deemed inappropriate by the principal investigator.

1.4.3 Testing Reagent

A. Evaluated Reagent

Reagent Name: Rapid SARS-CoV-2 Antigen Test Card

Packaging: 20 tests/kit

Storage: 2-30,°C sealed and stored in dry conditions

B. Comparator Method: RT-PCR method

Reagent Name: Detection Kit for 2019-nCoV (PCR-Fluorescence)

Manufacturer: Da An Gene Co., Ltd. of Sun Yat-sen University

Reagent Name: Diagnovital SARS-CoV-2 Real-Time PCR Kit

Manufacturer: A1 Life Sciences

Reagent Name: Allplex™ SARS-CoV-2 Assay

Manufacturer: Seegene Inc.

1.4.4 Quality Control Methods

1. Adopt unified technical methods in the same laboratory. Evaluation should be performed by the same testing personnel in strict accordance with the reagent operation procedures.
2. Indoor quality control should be carried out for each test.

1.4.5 Statistical Analysis Methods

1. The sensitivity, specificity, accuracy and confidence interval were calculated by comparing the detection results of Evaluated Reagent with RT-PCR (see Table 1).

Table 1 Comparison of test results for the Evaluated Reagent and RT-PCR

Evaluated Reagent Results	RT-PCR Results		Total
	Positive (+)	Negative (-)	
Positive (+)	a	b	a+b
Negative (-)	c	d	c+d
Total	a+c	b+d	a+b+c+d

Sensitivity (PPA) = $a/(a+c) * 100\%$;

Specificity (NPA) = $d/(b+d) * 100\%$;

Accuracy (OPA) = $(a+d)/(a+b+c+d) * 100\%$;

PPV = $a/(a+b) * 100\%$;

NPV = $d/(c+d) * 100\%$



MP Biomedicals Germany GmbH
Thueringer Straße 15
37269 Eschwege GERMANY
mpbio.com

2. Stratified statistics were made for positive samples with different onset days.

Table 2 Detection rate of positive samples with different onset days

	Number of Samples	Detection Rate
1-3 days		
4-7 days		

1.5 Clinical Trial Results and Analysis

1.5.1 Sample distribution statistics

A total of 1027 samples were used for the clinical trials, including 313 positive samples and 714 negative samples.

Table 3 Sample Distribution

	Nasopharyngeal Swab Samples	Nasal Swab Samples	Oropharyngeal Swab Samples	Total
Positive samples	102	105	106	313
Negative samples	464	125	125	714
Total	566	230	231	1027

1.5.2 Statistical analysis of datas

a. For nasopharyngeal swab samples

Table 4 Comparison of test results for the Evaluated Reagent and RT-PCR

Evaluated Reagent Results	RT-PCR Results		Total
	Positive (+)	Negative (-)	
Positive (+)	98	4	102
Negative (-)	4	460	464
Total	102	464	566

Sensitivity (PPA) = $98/102 \times 100\% = 96.08\%$ (95% CI: 92.31% - 99.85%)

Specificity (NPA) = $460/464 \times 100\% = 99.14\%$ (95% CI: 98.30% - 99.98%)

Accuracy (OPA) = $558/566 \times 100\% = 98.59\%$ (95% CI: 97.61% - 99.56%)

PPV = $98/102 \times 100\% = 96.08\%$ (95% CI: 92.31% - 99.85%)

NPV = $460/464 \times 100\% = 99.14\%$ (95% CI: 98.30% - 99.98%)

Sensitivity (PPA): $98/102 \times 100\% = 96.08\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9608 \pm 1.96 \times [0.9608 \times (1-0.9608) / 102]^{1/2}$
= [92.31%, 99.85%]

Specificity (NPA): $460/464 \times 100\% = 99.14\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9914 \pm 1.96 \times [0.9914 \times (1-0.9914) / 464]^{1/2}$
= [98.30%, 99.98%]



MP Biomedicals Germany GmbH
Thueringer Straße 15
37269 Eschwege GERMANY
mpbio.com

Accuracy (OPA): $558/566 \times 100\% = 98.59\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9859 \pm 1.96 \times [0.9859 \times (1-0.9859) / 566]^{1/2}$
= [97.61%, 99.56%]

PPV: $98/102 \times 100\% = 96.08\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9608 \pm 1.96 \times [0.9608 \times (1-0.9608) / 102]^{1/2}$
= [92.31%, 99.85%]

NPV: $460/464 \times 100\% = 99.14\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9914 \pm 1.96 \times [0.9914 \times (1-0.9914) / 464]^{1/2}$
= [98.30%, 99.98%]

b. For nasal swab samples

Table 5 Comparison of test results for the Evaluated Reagent and RT-PCR

Evaluated Reagent Results	RT-PCR Results		Total
	Positive (+)	Negative (-)	
Positive (+)	101	1	102
Negative (-)	4	124	128
Total	105	125	230

Sensitivity (PPA) = $101/105 \times 100\% = 96.19\%$ (95% CI: 92.53% - 99.85%)

Specificity (NPA) = $124/125 \times 100\% = 99.20\%$ (95% CI: 97.64% - 99.99%)

Accuracy (OPA) = $225/230 \times 100\% = 97.83\%$ (95% CI: 95.94% - 99.71%)

Positive Predictive Value (PPV) = $101/102 \times 100\% = 99.02\%$ (95% CI: 97.11% - 99.99%)

Negative Predictive Value (NPV) = $124/128 \times 100\% = 96.88\%$ (95% CI: 93.86% - 99.89%)

Sensitivity (PPA): $101/105 \times 100\% = 96.19\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9619 \pm 1.96 \times [0.9619 \times (1-0.9619) / 105]^{1/2}$
= [92.53%, 99.85%]

Specificity (NPA): $124/125 \times 100\% = 99.20\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9920 \pm 1.96 \times [0.9920 \times (1-0.9920) / 125]^{1/2}$
= [97.64%, 99.99%]

Accuracy (OPA): $225/230 \times 100\% = 97.83\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9783 \pm 1.96 \times [0.9783 \times (1-0.9783) / 230]^{1/2}$
= [95.94%, 99.71%]

PPV: $101/102 \times 100\% = 99.02\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9902 \pm 1.96 \times [0.9902 \times (1-0.9902) / 102]^{1/2}$
= [97.11%, 99.99%]

NPV: $124/128 \times 100\% = 96.88\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9688 \pm 1.96 \times [0.9688 \times (1-0.9688) / 128]^{1/2}$
= [93.86%, 99.89%]



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

c. For oropharyngeal swab samples

Table 6 Comparison of test results for the Evaluated Reagent and RT-PCR

Evaluated Reagent Results	RT-PCR Results		Total
	Positive (+)	Negative (-)	
Positive (+)	102	1	103
Negative (-)	4	124	128
Total	106	125	231

Sensitivity (PPA) = $102/106 \times 100\% = 96.23\%$ (95% CI: 92.60% - 99.85%)

Specificity (NPA) = $124/125 \times 100\% = 99.20\%$ (95% CI: 97.64% - 99.99%)

Accuracy (OPA) = $226/231 \times 100\% = 97.84\%$ (95% CI: 95.96% - 99.71%)

PPV = $102/103 \times 100\% = 99.03\%$ (95% CI: 97.14% - 99.99%)

NPV = $124/128 \times 100\% = 96.88\%$ (95% CI: 93.86% - 99.89%)

Sensitivity (PPA): $102/106 \times 100\% = 96.23\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9623 \pm 1.96 \times [0.9623 \times (1-0.9623) / 106]^{1/2}$
 $= [92.60\%, 99.85\%]$

Specificity (NPA): $124/125 \times 100\% = 99.20\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9920 \pm 1.96 \times [0.9920 \times (1-0.9920) / 125]^{1/2}$
 $= [97.64\%, 99.99\%]$

Accuracy (OPA): $226/231 \times 100\% = 97.84\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9784 \pm 1.96 \times [0.9784 \times (1-0.9784) / 231]^{1/2}$
 $= [95.96\%, 99.71\%]$

PPV: $102/103 \times 100\% = 99.03\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9903 \pm 1.96 \times [0.9903 \times (1-0.9903) / 103]^{1/2}$
 $= [97.14\%, 99.99\%]$

NPV: $124/128 \times 100\% = 96.88\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9688 \pm 1.96 \times [0.9688 \times (1-0.9688) / 128]^{1/2}$
 $= [93.86\%, 99.89\%]$



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

Table 7 Comparison of test results for the Evaluated Reagent and RT-PCR

Evaluated Reagent Results	RT-PCR Results		Total
	Positive (+)	Negative (-)	
Positive (+)	301	6	307
Negative (-)	12	708	720
Total	313	714	1027

Sensitivity (PPA) = $301/313 \times 100\% = 96.17\%$ (95% CI: 94.04% - 98.29%)

Specificity (NPA) = $708/714 \times 100\% = 99.16\%$ (95% CI: 98.49% - 99.83%)

Accuracy (OPA) = $1009/1027 \times 100\% = 98.25\%$ (95% CI: 97.44% - 99.05%)

PPV = $301/307 \times 100\% = 98.05\%$ (95% CI: 96.50% - 99.59%)

NPV = $708/720 \times 100\% = 98.33\%$ (95% CI: 97.40% - 99.27%)

Sensitivity (PPA): $301/313 \times 100\% = 96.17\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9617 \pm 1.96 \times [0.9617 \times (1-0.9617) / 313]^{1/2}$
 = [94.04%, 98.29%]

Specificity (NPA): $708/714 \times 100\% = 99.16\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9916 \pm 1.96 \times [0.9916 \times (1-0.9916) / 714]^{1/2}$
 = [98.49%, 99.83%]

Accuracy (OPA): $1009/1027 \times 100\% = 98.25\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9825 \pm 1.96 \times [0.9825 \times (1-0.9825) / 1027]^{1/2}$
 = [97.44%, 99.05%]

PPV: $301/307 \times 100\% = 98.05\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9805 \pm 1.96 \times [0.9805 \times (1-0.9805) / 307]^{1/2}$
 = [96.50%, 99.59%]

NPV: $708/720 \times 100\% = 98.33\%$

95% CI: $p \pm 1.96 \times [p(1-p) / n]^{1/2} = 0.9833 \pm 1.96 \times [0.9833 \times (1-0.9833) / 720]^{1/2}$
 = [97.40%, 99.27%]

e. Detection rate of positive samples with different onset days

Table 8 Detection rate of positive samples with different onset days

	Number of Samples	Detection Rate
1-3 days	165	161/165, 97.58%
4-7 days	148	140/148, 94.59%

Detection rate of 1-3 days = $161/165 \times 100\% = 97.58\%$ (95% CI: 95.23% - 99.92%)

Detection rate of 4-7 days = $140/148 \times 100\% = 94.59\%$ (95% CI: 90.95% - 98.24%)



MP Biomedicals Germany GmbH
Thueringer Straße 15
37269 Eschwege GERMANY
mpbio.com

2 Appendix

Case Report Forms

CONFIDENTIAL



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

Appendix
Case Report Form

No.	Patient Number	Sex	Age	Sample Type	Clinical Diagnosis Background Info	Evaluated Reagent Results	Nucleic Acid Test Results	Nucleic Acid Testing Value (Ct)	Disease Onset Time
1	00015164	F	42	Nasopharyngeal Swab	COVID-19	+	+	25	5 days
2	00653145	F	65	Nasopharyngeal Swab	COVID-19	+	+	29	1 day
3	00024862	F	20	Nasopharyngeal Swab	COVID-19	+	+	21.5	3 days
4	00041020	F	38	Nasopharyngeal Swab	COVID-19	+	+	21	1 day
5	00079066	F	72	Nasopharyngeal Swab	COVID-19	+	+	22	3 days
6	00085102	M	66	Nasopharyngeal Swab	COVID-19	+	+	24.3	2 days
7	00070960	F	37	Nasopharyngeal Swab	COVID-19	+	+	23.3	2 days
8	00014950	M	38	Nasopharyngeal Swab	COVID-19	-	+	27	4 days
9	00064687	F	20	Nasopharyngeal Swab	COVID-19	+	+	22.6	2 days
10	00045019	F	26	Nasopharyngeal Swab	COVID-19	+	+	28	7 days
11	00034165	F	53	Nasopharyngeal Swab	COVID-19	+	+	23	1 day
12	00045147	F	86	Nasopharyngeal Swab	COVID-19	+	+	21	3 days
13	00065219	M	56	Nasopharyngeal Swab	COVID-19	+	+	22	4 days
14	00049469	M	58	Nasopharyngeal Swab	COVID-19	+	+	24.2	2 days
15	00085476	M	10	Nasopharyngeal Swab	COVID-19	+	+	23.1	3 days
16	00055367	M	34	Nasopharyngeal Swab	COVID-19	+	+	28	1 day
17	00065069	F	39	Nasopharyngeal Swab	COVID-19	+	+	27.4	2 days
18	00090180	F	47	Nasopharyngeal Swab	COVID-19	+	+	22	2 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

19	00012761	F	61	Nasopharyngeal Swab	COVID-19	+	+	27	6 days
20	00056782	M	30	Nasopharyngeal Swab	COVID-19	+	+	25.5	4 days
21	00043045	F	15	Nasopharyngeal Swab	COVID-19	+	+	23.6	1 day
22	00055198	F	47	Nasopharyngeal Swab	COVID-19	+	+	25.2	2 days
23	00032371	F	29	Nasopharyngeal Swab	COVID-19	+	+	26	2 days

24	00099760	F	49	Nasopharyngeal Swab	COVID-19	+	+	27.1	2 days
25	00056207	F	43	Nasopharyngeal Swab	COVID-19	+	+	31.2	4 days
26	00042236	M	58	Nasopharyngeal Swab	COVID-19	+	+	29	6 days
27	00089182	F	18	Nasopharyngeal Swab	COVID-19	+	+	24	2 days
28	00052691	F	13	Nasopharyngeal Swab	COVID-19	+	+	25	2 days
29	00036306	F	50	Nasopharyngeal Swab	COVID-19	+	+	22.5	2 days
30	00058181	F	55	Nasopharyngeal Swab	COVID-19	+	+	24.4	1 day
31	00024549	M	48	Nasopharyngeal Swab	COVID-19	-	+	26.5	3 days
32	00082479	F	74	Nasopharyngeal Swab	COVID-19	+	+	26.3	4 days
33	20201026003	F	36	Nasopharyngeal Swab	COVID-19	+	+	28.1	2 days
34	20201026006	M	25	Nasopharyngeal Swab	COVID-19	+	+	29.0	4 days
35	20201026010	F	34	Nasopharyngeal Swab	COVID-19	+	+	28.9	1 days
36	20201026013	M	19	Nasopharyngeal Swab	COVID-19	+	+	26.7	2 days
37	20201026018	M	28	Nasopharyngeal Swab	COVID-19	+	+	32.8	4 days
38	20201026019	M	45	Nasopharyngeal Swab	COVID-19	+	+	34.0	3 days
39	20201026020	F	34	Nasopharyngeal Swab	COVID-19	+	+	26.7	2 days
40	20201026021	M	55	Nasopharyngeal Swab	COVID-19	+	+	32.8	5 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

41	20201026022	M	32	Nasopharyngeal Swab	COVID-19	+	+	34.0	3 days
42	20201026027	F	22	Nasopharyngeal Swab	COVID-19	+	+	26.7	5 days
43	20201026034	F	28	Nasopharyngeal Swab	COVID-19	+	+	31.0	3 days
44	20201026037	F	35	Nasopharyngeal Swab	COVID-19	+	+	29.3	4 days
45	20201026043	M	56	Nasopharyngeal Swab	COVID-19	+	+	25.3	1 days
46	20201026048	M	61	Nasopharyngeal Swab	COVID-19	+	+	25.0	5 days
47	20201026123	F	23	Nasopharyngeal Swab	COVID-19	+	+	24.87	2 days
48	20201026126	F	21	Nasopharyngeal Swab	COVID-19	+	+	29.49	1 day
49	20201026134	F	34	Nasopharyngeal Swab	COVID-19	+	+	27.42	2 days
50	20201026159	M	27	Nasopharyngeal Swab	COVID-19	+	+	31.25	3 days
51	20201026057	M	26	Nasopharyngeal Swab	COVID-19	+	+	25.6	3 days
52	20201026063	F	56	Nasopharyngeal Swab	COVID-19	+	+	26.5	2 days

53	20201026068	F	32	Nasopharyngeal Swab	COVID-19	+	+	26.8	5 days
54	20201026072	F	27	Nasopharyngeal Swab	COVID-19	+	+	23.6	3 days
55	20201026077	M	43	Nasopharyngeal Swab	COVID-19	+	+	25.0	1 days
56	20201026078	M	36	Nasopharyngeal Swab	COVID-19	+	+	26.8	2 days
57	20201026088	F	43	Nasopharyngeal Swab	COVID-19	+	+	23.6	5 days
58	20201026079	F	34	Nasopharyngeal Swab	COVID-19	+	+	25.0	6 days
59	20201026080	M	56	Nasopharyngeal Swab	COVID-19	+	+	26.8	4 days
60	20201026082	M	54	Nasopharyngeal Swab	COVID-19	+	+	20.7	2 days
61	20201026085	M	47	Nasopharyngeal Swab	COVID-19	+	+	26.1	3 days
62	20201026093	F	46	Nasopharyngeal Swab	COVID-19	+	+	20.0	5 days
63	20201026096	M	22	Nasopharyngeal Swab	COVID-19	+	+	27.1	6 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

64	20201026097	F	32	Nasopharyngeal Swab	COVID-19	+	+	22.8	2 days
65	20201026128	M	48	Nasopharyngeal Swab	COVID-19	+	+	27.43	3 days
66	20201026129	F	43	Nasopharyngeal Swab	COVID-19	+	+	29.79	1 days
67	20201026132	F	19	Nasopharyngeal Swab	COVID-19	+	+	30.58	1 days
68	20201026138	F	22	Nasopharyngeal Swab	COVID-19	+	+	30.26	2 days
69	201202001	M	76	Nasopharyngeal Swab	COVID-19	+	+	20.16	5 days
70	201202002	M	75	Nasopharyngeal Swab	COVID-19	+	+	27.78	4 days
71	201202003	F	82	Nasopharyngeal Swab	COVID-19	+	+	22.73	4 days
72	201202004	F	66	Nasopharyngeal Swab	COVID-19	+	+	17.6	3 days
73	201202005	M	56	Nasopharyngeal Swab	COVID-19	+	+	27.21	7 days
74	201202006	F	59	Nasopharyngeal Swab	COVID-19	+	+	20.24	2 days
75	201202007	F	68	Nasopharyngeal Swab	COVID-19	+	+	18.58	1 day
76	201202008	F	62	Nasopharyngeal Swab	COVID-19	+	+	23.78	4 days
77	201202009	M	83	Nasopharyngeal Swab	COVID-19	+	+	24.06	5 days
78	201202010	M	52	Nasopharyngeal Swab	COVID-19	+	+	16.04	3 days
79	201202011	M	73	Nasopharyngeal Swab	COVID-19	-	+	25.06	3 days
80	201202012	F	24	Nasopharyngeal Swab	COVID-19	+	+	17.84	3 days
81	201202013	F	59	Nasopharyngeal Swab	COVID-19	+	+	29.5	7 days
82	201202014	F	66	Nasopharyngeal Swab	COVID-19	+	+	33.92	5 days
83	201202015	F	36	Nasopharyngeal Swab	COVID-19	+	+	26.98	4 days
84	201202016	F	24	Nasopharyngeal Swab	COVID-19	+	+	33.14	3 days
85	201202017	F	27	Nasopharyngeal Swab	COVID-19	+	+	17.94	6 days
86	201202018	M	27	Nasopharyngeal Swab	COVID-19	+	+	20.28	2 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

87	201202019	M	23	Nasopharyngeal Swab	COVID-19	+	+	20.4	4 days
88	201202020	F	64	Nasopharyngeal Swab	COVID-19	+	+	16.64	3 days
89	201202021	M	67	Nasopharyngeal Swab	COVID-19	+	+	21.76	6 days
90	201202022	F	32	Nasopharyngeal Swab	COVID-19	+	+	22.86	1 day
91	201202023	M	79	Nasopharyngeal Swab	COVID-19	-	+	27.41	5 days
92	201202024	F	59	Nasopharyngeal Swab	COVID-19	+	+	22.17	4 days
93	201202025	M	64	Nasopharyngeal Swab	COVID-19	+	+	24.45	5 days
94	201202026	F	62	Nasopharyngeal Swab	COVID-19	+	+	25.77	7 days
95	201202027	F	63	Nasopharyngeal Swab	COVID-19	+	+	24.84	5 days
96	201202028	F	55	Nasopharyngeal Swab	COVID-19	+	+	24.25	2 days
97	201202029	F	58	Nasopharyngeal Swab	COVID-19	+	+	27.86	4 days
98	201202030	M	57	Nasopharyngeal Swab	COVID-19	+	+	15.98	7 days
99	201202031	F	57	Nasopharyngeal Swab	COVID-19	+	+	31.35	3 days
100	201202032	F	60	Nasopharyngeal Swab	COVID-19	+	+	30.1	4 days
101	201202033	F	30	Nasopharyngeal Swab	COVID-19	+	+	33.6	6 days
102	201202034	F	66	Nasopharyngeal Swab	COVID-19	+	+	31.33	5 days
103	2020037819	F	22	Nasal Swab	COVID-19	+	+	28.98	7 days
104	2020037657	F	37	Nasal Swab	COVID-19	+	+	20.21	3 days
105	2020037480	M	33	Nasal Swab	COVID-19	+	+	26.94	1 day
106	2020037273	M	17	Nasal Swab	COVID-19	+	+	26.10	6 days
107	2020036912	M	47	Nasal Swab	COVID-19	+	+	22.71	7 days
108	2020036814	F	20	Nasal Swab	COVID-19	+	+	28.42	2 days
109	2020036496	M	16	Nasal Swab	COVID-19	+	+	22.50	3 days



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

110	2020036461	M	16	Nasal Swab	COVID-19	+	+	24.00	4 days
111	2020035999	M	41	Nasal Swab	COVID-19	+	+	27.07	1 day
112	2020035992	M	29	Nasal Swab	COVID-19	+	+	27.54	2 days
113	2020035758	F	41	Nasal Swab	COVID-19	+	+	24.81	1 day
114	2020035710	F	50	Nasal Swab	COVID-19	+	+	28.74	7 days
115	2020035661	M	23	Nasal Swab	COVID-19	+	+	20.85	3 days
116	2020035648	M	37	Nasal Swab	COVID-19	+	+	22.97	5 days
117	2020035443	F	62	Nasal Swab	COVID-19	+	+	28.37	4 days
118	2020035164	F	23	Nasal Swab	COVID-19	-	+	29.48	6 days
119	2020034977	M	30	Nasal Swab	COVID-19	+	+	20.80	7 days
120	2020034774	F	40	Nasal Swab	COVID-19	+	+	27.40	4 days
121	2020034711	M	18	Nasal Swab	COVID-19	+	+	23.05	5 days
122	2020034169	M	32	Nasal Swab	COVID-19	+	+	23.23	3 days
123	2020033951	M	64	Nasal Swab	COVID-19	+	+	21.42	4 days
124	2020033575	F	27	Nasal Swab	COVID-19	+	+	30.65	5 days
125	2020033077	F	43	Nasal Swab	COVID-19	+	+	26.77	1 day
126	2020032904	M	31	Nasal Swab	COVID-19	+	+	30.87	6 days
127	2020032594	M	54	Nasal Swab	COVID-19	+	+	20.87	3 days
128	2020032584	F	48	Nasal Swab	COVID-19	+	+	30.74	2 days
129	2020032249	M	24	Nasal Swab	COVID-19	+	+	23.82	5 days
130	2020032204	M	65	Nasal Swab	COVID-19	+	+	24.89	3 days
131	2020032015	M	50	Nasal Swab	COVID-19	+	+	21.58	6 days
132	2020031510	M	35	Nasal Swab	COVID-19	+	+	25.72	5 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

133	2020031306	M	27	Nasal Swab	COVID-19	+	+	26.18	3 days
134	210085230	F	57	Nasal Swab	COVID-19	+	+	33.18	5 days
135	210085383	M	47	Nasal Swab	COVID-19	+	+	25.02	6 days
136	210085421	F	42	Nasal Swab	COVID-19	+	+	19.49	1 day
137	210085431	F	44	Nasal Swab	COVID-19	+	+	27.92	6 days
138	210085435	F	40	Nasal Swab	COVID-19	+	+	28.80	7 days
139	210085512	M	63	Nasal Swab	COVID-19	+	+	31.50	3 days

140	210085544	M	43	Nasal Swab	COVID-19	+	+	30.54	7 days
141	210085560	F	31	Nasal Swab	COVID-19	+	+	32.16	1 day
142	210085593	M	46	Nasal Swab	COVID-19	+	+	32.07	7 days
143	210085606	M	36	Nasal Swab	COVID-19	+	+	31.39	4 days
144	210085609	F	57	Nasal Swab	COVID-19	+	+	23.52	3 days
145	210085617	M	42	Nasal Swab	COVID-19	+	+	26.13	2 days
146	210085657	F	63	Nasal Swab	COVID-19	+	+	34.77	4 days
147	210085676	F	45	Nasal Swab	COVID-19	+	+	33.29	4 days
148	210085725	F	24	Nasal Swab	COVID-19	+	+	18.57	7 days
149	210085733	M	67	Nasal Swab	COVID-19	+	+	18.62	7 days
150	210085775	M	29	Nasal Swab	COVID-19	+	+	23.19	1 day
151	210085821	F	53	Nasal Swab	COVID-19	+	+	23.76	2 days
152	210085834	M	49	Nasal Swab	COVID-19	+	+	25.22	1 day
153	210085838	F	59	Nasal Swab	COVID-19	+	+	25.86	6 days
154	210085869	M	55	Nasal Swab	COVID-19	+	+	29.95	5 days
155	210085978	M	36	Nasal Swab	COVID-19	+	+	28.82	4 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

156	210086041	F	65	Nasal Swab	COVID-19	+	+	28.87	1 day
157	210086048	M	23	Nasal Swab	COVID-19	+	+	21.52	1 day
158	210086049	F	33	Nasal Swab	COVID-19	+	+	23.43	6 days
159	210086051	M	27	Nasal Swab	COVID-19	+	+	33.80	3 days
160	210086161	M	57	Nasal Swab	COVID-19	+	+	33.47	5 days
161	210086167	F	45	Nasal Swab	COVID-19	+	+	34.07	1 day
162	210086172	F	38	Nasal Swab	COVID-19	+	+	27.10	3 days
163	210086303	F	37	Nasal Swab	COVID-19	+	+	23.78	5 days
164	210086400	F	63	Nasal Swab	COVID-19	-	+	34.70	7 days
165	210086423	F	61	Nasal Swab	COVID-19	+	+	32.27	4 days
166	210086462	M	49	Nasal Swab	COVID-19	+	+	31.97	4 days
167	210086525	M	43	Nasal Swab	COVID-19	+	+	32.75	6 days
168	210086555	F	38	Nasal Swab	COVID-19	+	+	33.18	7 days

169	210086559	M	40	Nasal Swab	COVID-19	+	+	32.93	3 days
170	210086622	M	66	Nasal Swab	COVID-19	+	+	31.77	4 days
171	210086647	F	59	Nasal Swab	COVID-19	+	+	31.25	5 days
172	210086662	F	60	Nasal Swab	COVID-19	+	+	30.53	3 days
173	210086703	F	39	Nasal Swab	COVID-19	+	+	31.57	2 days
174	210086712	M	57	Nasal Swab	COVID-19	+	+	32.94	3 days
175	210086739	M	67	Nasal Swab	COVID-19	+	+	28.83	6 days
176	210086762	F	66	Nasal Swab	COVID-19	+	+	22.75	2 days
177	210086792	M	24	Nasal Swab	COVID-19	+	+	19.01	3 days
178	210086825	F	35	Nasal Swab	COVID-19	+	+	24.25	7 days



MP Biomedicals Germany GmbH
 Thüringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

179	210086858	M	29	Nasal Swab	COVID-19	+	+	27.56	7 days
180	210086896	F	27	Nasal Swab	COVID-19	+	+	25.73	4 days
181	210086902	M	47	Nasal Swab	COVID-19	+	+	31.37	4 days
182	210086908	M	66	Nasal Swab	COVID-19	+	+	31.93	3 days
183	210086918	F	56	Nasal Swab	COVID-19	+	+	30.07	2 days
184	210087237	M	42	Nasal Swab	COVID-19	+	+	26.61	3 days
185	210087242	F	33	Nasal Swab	COVID-19	+	+	32.27	6 days
186	210087256	M	59	Nasal Swab	COVID-19	+	+	27.53	3 days
187	210087319	M	29	Nasal Swab	COVID-19	+	+	30.43	2 days
188	210087397	F	19	Nasal Swab	COVID-19	+	+	22.36	1 day
189	210087399	M	35	Nasal Swab	COVID-19	+	+	29.47	5 days
190	210087417	F	28	Nasal Swab	COVID-19	+	+	25.65	2 days
191	210087422	F	23	Nasal Swab	COVID-19	+	+	28.06	4 days
192	210087427	F	63	Nasal Swab	COVID-19	+	+	31.03	3 days
193	210087431	M	54	Nasal Swab	COVID-19	+	+	23.15	6 days
194	210087453	M	38	Nasal Swab	COVID-19	-	+	34.64	7 days
195	210087461	F	40	Nasal Swab	COVID-19	+	+	28.49	2 days
196	210087482	F	32	Nasal Swab	COVID-19	+	+	23.62	4 days
197	210087487	M	25	Nasal Swab	COVID-19	+	+	26.25	3 days
198	210087493	F	41	Nasal Swab	COVID-19	+	+	29.91	4 days
199	210087511	F	64	Nasal Swab	COVID-19	+	+	30.67	2 days
200	210087534	M	47	Nasal Swab	COVID-19	-	+	33.46	3 days
201	210087540	F	33	Nasal Swab	COVID-19	+	+	31.24	5 days



MP Biomedicals Germany GmbH
 Thüringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

202	210087549	M	52	Nasal Swab	COVID-19	+	+	20.73	3 days
203	210087553	F	26	Nasal Swab	COVID-19	+	+	31.48	2 days
204	210087613	M	59	Nasal Swab	COVID-19	+	+	25.93	7 days
205	210087637	F	34	Nasal Swab	COVID-19	+	+	24.58	3 days
206	210087656	F	61	Nasal Swab	COVID-19	+	+	21.84	2 days
207	210087854	M	44	Nasal Swab	COVID-19	+	+	22.71	1 day
208	2020037203	M	63	Oropharyngeal Swab	COVID-19	+	+	23.63	5 days
209	2020037143	M	22	Oropharyngeal Swab	COVID-19	+	+	26.39	1 day
210	2020036948	F	24	Oropharyngeal Swab	COVID-19	+	+	22.99	7 days
211	2020036403	F	64	Oropharyngeal Swab	COVID-19	+	+	20.92	3 days
212	2020036386	M	58	Oropharyngeal Swab	COVID-19	+	+	24.19	2 days
213	2020036220	M	26	Oropharyngeal Swab	COVID-19	+	+	31.69	7 days
214	2020036135	M	58	Oropharyngeal Swab	COVID-19	+	+	30.36	4 days
215	2020036127	F	25	Oropharyngeal Swab	COVID-19	+	+	29.54	2 days
216	2020035970	F	58	Oropharyngeal Swab	COVID-19	+	+	21.63	1 day
217	2020035395	M	15	Oropharyngeal Swab	COVID-19	+	+	24.31	4 days
218	2020035239	M	38	Oropharyngeal Swab	COVID-19	+	+	28.96	7 days
219	2020035081	M	43	Oropharyngeal Swab	COVID-19	+	+	30.95	2 days
220	2020034843	M	52	Oropharyngeal Swab	COVID-19	+	+	31.63	5 days
221	2020034648	F	23	Oropharyngeal Swab	COVID-19	+	+	25.58	3 days
222	2020034553	F	42	Oropharyngeal Swab	COVID-19	+	+	30.47	4 days
223	2020034488	M	25	Oropharyngeal Swab	COVID-19	+	+	29.31	6 days
224	2020034171	F	31	Oropharyngeal Swab	COVID-19	+	+	29.13	1 day



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

225	2020034049	M	51	Oropharyngeal Swab	COVID-19	+	+	28.34	5 days
226	2020033543	F	30	Oropharyngeal Swab	COVID-19	+	+	24.03	3 days

227	2020033110	M	37	Oropharyngeal Swab	COVID-19	+	+	21.88	5 days
228	2020032950	M	30	Oropharyngeal Swab	COVID-19	+	+	25.00	4 days
229	2020032615	F	34	Oropharyngeal Swab	COVID-19	+	+	24.79	4 days
230	2020032503	F	38	Oropharyngeal Swab	COVID-19	+	+	30.39	7 days
231	2020032293	M	50	Oropharyngeal Swab	COVID-19	+	+	21.04	3 days
232	2020032229	F	49	Oropharyngeal Swab	COVID-19	+	+	28.08	6 days
233	2020032051	M	49	Oropharyngeal Swab	COVID-19	+	+	21.83	3 days
234	2020032040	F	53	Oropharyngeal Swab	COVID-19	+	+	24.89	1 day
235	2020031733	M	16	Oropharyngeal Swab	COVID-19	-	+	31.39	6 days
236	2020031578	M	58	Oropharyngeal Swab	COVID-19	+	+	25.31	5 days
237	2020031544	F	31	Oropharyngeal Swab	COVID-19	+	+	28.92	2 days
238	2020031480	F	55	Oropharyngeal Swab	COVID-19	+	+	29.66	6 days
239	210087857	M	32	Oropharyngeal Swab	COVID-19	+	+	25.12	3 days
240	210087860	F	39	Oropharyngeal Swab	COVID-19	+	+	21.54	3 days
241	210087862	F	52	Oropharyngeal Swab	COVID-19	+	+	21.03	3 days
242	210087864	M	32	Oropharyngeal Swab	COVID-19	+	+	22.86	4 days
243	210087866	F	46	Oropharyngeal Swab	COVID-19	+	+	24.33	1 day
244	210087867	F	67	Oropharyngeal Swab	COVID-19	+	+	33.16	5 days
245	210087868	M	34	Oropharyngeal Swab	COVID-19	+	+	27.19	3 days
246	210087870	F	25	Oropharyngeal Swab	COVID-19	+	+	22.48	1 day
247	210087873	F	56	Oropharyngeal Swab	COVID-19	+	+	28.47	2 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

248	210087874	M	41	Oropharyngeal Swab	COVID-19	+	+	23.69	5 days
249	210087876	F	65	Oropharyngeal Swab	COVID-19	+	+	31.18	4 days
250	210087878	M	38	Oropharyngeal Swab	COVID-19	+	+	22.49	2 days
251	210087881	F	44	Oropharyngeal Swab	COVID-19	+	+	24.72	3 days
252	210087884	M	67	Oropharyngeal Swab	COVID-19	+	+	33.51	5 days
253	210087887	F	24	Oropharyngeal Swab	COVID-19	-	+	28.46	3 days
254	210087891	F	27	Oropharyngeal Swab	COVID-19	+	+	27.19	3 days
255	210087894	F	57	Oropharyngeal Swab	COVID-19	+	+	22.58	5 days

256	210087895	M	63	Oropharyngeal Swab	COVID-19	+	+	27.97	3 days
257	210087897	M	40	Oropharyngeal Swab	COVID-19	+	+	25.54	1 day
258	210087902	F	32	Oropharyngeal Swab	COVID-19	+	+	23.62	2 days
259	210087904	F	37	Oropharyngeal Swab	COVID-19	+	+	32.16	4 days
260	210087907	M	45	Oropharyngeal Swab	COVID-19	+	+	26.83	7 days
261	210087908	F	61	Oropharyngeal Swab	COVID-19	+	+	27.15	3 days
262	210087910	M	31	Oropharyngeal Swab	COVID-19	+	+	31.24	4 days
263	210087914	F	38	Oropharyngeal Swab	COVID-19	+	+	29.95	3 days
264	210087916	F	63	Oropharyngeal Swab	COVID-19	+	+	24.16	3 days
265	210087918	M	21	Oropharyngeal Swab	COVID-19	+	+	25.83	3 days
266	210087921	F	46	Oropharyngeal Swab	COVID-19	+	+	22.52	1 day
267	210087923	F	23	Oropharyngeal Swab	COVID-19	+	+	24.49	6 days
268	210087925	M	35	Oropharyngeal Swab	COVID-19	+	+	26.57	5 days
269	210087927	F	69	Oropharyngeal Swab	COVID-19	-	+	34.33	7 days
270	210087930	F	55	Oropharyngeal Swab	COVID-19	+	+	28.14	2 days



MP Biomedicals Germany GmbH
 Thüringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

271	210087932	F	42	Oropharyngeal Swab	COVID-19	+	+	29.53	3 days
272	210087935	M	30	Oropharyngeal Swab	COVID-19	+	+	28.95	1 day
273	210087937	F	49	Oropharyngeal Swab	COVID-19	+	+	26.72	5 days
274	210087938	M	46	Oropharyngeal Swab	COVID-19	+	+	32.08	3 days
275	210087941	F	28	Oropharyngeal Swab	COVID-19	+	+	31.15	4 days
276	210087943	M	41	Oropharyngeal Swab	COVID-19	+	+	26.74	2 days
277	210087948	M	56	Oropharyngeal Swab	COVID-19	+	+	22.81	5 days
278	210087950	F	37	Oropharyngeal Swab	COVID-19	+	+	30.47	7 days
279	210087952	F	24	Oropharyngeal Swab	COVID-19	+	+	23.24	3 days
280	210087954	M	43	Oropharyngeal Swab	COVID-19	-	+	31.24	6 days
281	210087955	M	52	Oropharyngeal Swab	COVID-19	+	+	23.35	5 days
282	210087957	M	18	Oropharyngeal Swab	COVID-19	+	+	25.13	1 day
283	210087959	M	44	Oropharyngeal Swab	COVID-19	+	+	25.86	5 days
284	210087962	F	37	Oropharyngeal Swab	COVID-19	+	+	24.87	3 days

285	210087965	F	28	Oropharyngeal Swab	COVID-19	+	+	29.49	2 days
286	210087967	M	47	Oropharyngeal Swab	COVID-19	+	+	27.42	4 days
287	210087970	F	58	Oropharyngeal Swab	COVID-19	+	+	31.25	3 days
288	210087972	F	62	Oropharyngeal Swab	COVID-19	+	+	25.06	2 days
289	210087975	M	40	Oropharyngeal Swab	COVID-19	+	+	26.43	5 days
290	210087983	F	29	Oropharyngeal Swab	COVID-19	+	+	26.78	3 days
291	210087986	M	54	Oropharyngeal Swab	COVID-19	+	+	23.36	6 days
292	210087989	F	33	Oropharyngeal Swab	COVID-19	+	+	25.19	4 days
293	210087992	M	37	Oropharyngeal Swab	COVID-19	+	+	22.18	2 days



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

294	210087995	M	42	Oropharyngeal Swab	COVID-19	+	+	23.06	5 days
295	210087998	M	50	Oropharyngeal Swab	COVID-19	+	+	25.28	3 days
296	210088005	F	72	Oropharyngeal Swab	COVID-19	+	+	26.15	4 days
297	210088006	F	53	Oropharyngeal Swab	COVID-19	+	+	20.87	3 days
298	210088013	M	46	Oropharyngeal Swab	COVID-19	+	+	29.35	6 days
299	210088016	F	24	Oropharyngeal Swab	COVID-19	+	+	25.21	5 days
300	210088019	M	59	Oropharyngeal Swab	COVID-19	+	+	20.84	2 day
301	210088022	F	28	Oropharyngeal Swab	COVID-19	+	+	27.19	3 days
302	210088025	F	22	Oropharyngeal Swab	COVID-19	+	+	29.83	4 days
303	210088028	M	34	Oropharyngeal Swab	COVID-19	+	+	19.98	6 days
304	210088031	F	46	Oropharyngeal Swab	COVID-19	+	+	20.33	2 days
305	210088033	M	39	Oropharyngeal Swab	COVID-19	+	+	23.04	3 days
306	210088036	F	45	Oropharyngeal Swab	COVID-19	+	+	26.31	1 day
307	210088039	M	24	Oropharyngeal Swab	COVID-19	+	+	29.46	3 days
308	210088042	F	60	Oropharyngeal Swab	COVID-19	+	+	20.19	3 days
309	210088045	M	45	Oropharyngeal Swab	COVID-19	+	+	20.33	5 days
310	210088047	F	64	Oropharyngeal Swab	COVID-19	+	+	22.79	3 days
311	210088052	M	26	Oropharyngeal Swab	COVID-19	+	+	31.36	4 days
312	210088058	F	44	Oropharyngeal Swab	COVID-19	+	+	24.34	6 days
313	210088071	F	32	Oropharyngeal Swab	COVID-19	+	+	27.85	3 days
314	00038365	M	70	Nasopharyngeal Swab	AIDS	-	-	> 40	/
315	00076707	M	39	Nasopharyngeal Swab	AIDS	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

316	00040634	M	51	Nasopharyngeal Swab	AIDS	-	-	> 40	/
317	00039635	F	38	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
318	QY005248	M	65	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
319	00085216	M	44	Nasopharyngeal Swab	Tuberculous Meningitis	-	-	> 40	/
320	00048479	M	22	Nasopharyngeal Swab	Cough to be tested	-	-	> 40	/
321	00092844	M	34	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
322	00034321	M	25	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
323	00043099	M	25	Nasopharyngeal Swab	Tuberculosis re-treatment	-	-	> 40	/
324	00062236	M	56	Nasopharyngeal Swab	Tuberculosis, Culture-positive, Tracheal Tuberculosis	+	-	> 40	/
325	00012535	F	31	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
326	00060272	F	13	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
327	00039324	M	54	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
328	00025899	M	74	Nasopharyngeal Swab	Fever to be tested, Obsolete Tuberculosis	-	-	> 40	/
329	00654962	M	75	Nasopharyngeal Swab	Fever to be tested, Hand Furuncle with Infection?	-	-	> 40	/
330	00834102	F	37	Nasopharyngeal Swab	Fever to be tested: 1) bacterial? 2) viral?, Hypothyroidism	+	-	> 40	/
331	00454634	F	20	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/
332	00619267	F	69	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

333	00126784	F	30	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/
334	00449427	M	26	Nasopharyngeal Swab	Fever to be tested: Lung Infection?	-	-	> 40	/

335	00383918	F	38	Nasopharyngeal Swab	Pneumonia recovery period	-	-	> 40	/
336	00964723	M	9	Nasopharyngeal Swab	URTI?	-	-	> 40	/
337	00306115	F	49	Nasopharyngeal Swab	URTI?	-	-	> 40	/
338	00130338	M	30	Nasopharyngeal Swab	URTI?	-	-	> 40	/
339	00411660	M	37	Nasopharyngeal Swab	URTI?	-	-	> 40	/
340	00778439	F	81	Nasopharyngeal Swab	Trochanteric fracture of left femur	-	-	> 40	/
341	00583893	M	72	Nasopharyngeal Swab	COPD with acute exacerbation	-	-	> 40	/
342	00454385	F	29	Nasopharyngeal Swab	Early Intrauterine Pregnancy	-	-	> 40	/
343	00974308	F	53	Nasopharyngeal Swab	Community Acquired Pneumonia	-	-	> 40	/
344	00429842	F	21	Nasopharyngeal Swab	Mycoplasma Pneumoniae Pneumonia, Community Acquired Pneumonia	-	-	> 40	/
345	00231585	M	38	Nasopharyngeal Swab	Community Acquired Pneumonia, non-severe	-	-	> 40	/
346	00414097	M	66	Nasopharyngeal Swab	COPD with acute exacerbation	-	-	> 40	/
347	00505685	M	60	Nasopharyngeal Swab	Pulmonary Heart Disease	+	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

348	00514175	F	82	Nasopharyngeal Swab	Fell and hurt	-	-	> 40	/
349	00534501	F	39	Nasopharyngeal Swab	Early Pregnancy Miscarriage	-	-	> 40	/
350	00489143	F	48	Nasopharyngeal Swab	AIDS	-	-	> 40	/
351	00765061	M	64	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
352	00527015	M	38	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
353	00567376	M	33	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
354	00850446	M	80	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
355	00844574	M	73	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
356	00724394	M	25	Nasopharyngeal Swab	AIDS	-	-	> 40	/

357	00323014	F	54	Nasopharyngeal Swab	AIDS	-	-	> 40	/
358	00527445	M	77	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
359	00973853	M	60	Nasopharyngeal Swab	AIDS	-	-	> 40	/
360	00489735	M	25	Nasopharyngeal Swab	AIDS	-	-	> 40	/
361	00883964	F	14	Nasopharyngeal Swab	Cervical and abdominal lymph node TB	-	-	> 40	/
362	00393257	M	78	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
363	00665674	M	84	Nasopharyngeal Swab	Tuberculosis, post-operative multiple-drug resistance	-	-	> 40	/
364	00973975	M	57	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
365	00369045	M	65	Nasopharyngeal Swab	Chronic HBV	-	-	> 40	/
366	00081637	M	44	Nasopharyngeal Swab	AIDS, Tuberculosis, Drug-induced liver damage	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

367	00013753	M	37	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
368	00098670	F	28	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
369	00043616	M	21	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
370	00057174	M	40	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
371	00034718	M	75	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
372	00064734	M	35	Nasopharyngeal Swab	Abdominal distention to be tested	-	-	> 40	/
373	00081698	F	52	Nasopharyngeal Swab	Liver Cirrhosis	-	-	> 40	/
374	00023776	M	49	Nasopharyngeal Swab	HBV	-	-	> 40	/
375	00049609	M	51	Nasopharyngeal Swab	Upper Gastrointestinal Bleeding	-	-	> 40	/
376	00060356	F	63	Nasopharyngeal Swab	Paravertebral Abscess	-	-	> 40	/
377	00056728	M	63	Nasopharyngeal Swab	AIDS, Lung Infection, Headache to be tested	-	-	> 40	/
378	00026576	M	23	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
379	00032806	M	25	Nasopharyngeal Swab	Cough to be tested	-	-	> 40	/
380	00036913	F	30	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
381	00020793	M	66	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thüringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

382	00076534	F	31	Nasopharyngeal Swab	Secondary Tuberculosis (re-treatment, multiple drug resistance) smear- and culture-positive, Bronchial Tuberculosis, Drug-Induced Liver Injury, Vitamin B6 deficiency	-	-	> 40	/
383	00051375	M	55	Nasopharyngeal Swab	AIDS, fever to be tested	-	-	> 40	/
384	00052681	M	35	Nasopharyngeal Swab	AIDS, fever to be tested	-	-	> 40	/
385	00013626	F	23	Nasopharyngeal Swab	Fever	-	-	> 40	/
386	00046289	M	50	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
387	00059952	M	39	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
388	00019759	M	40	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
389	00038800	M	81	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
390	00055841	M	64	Nasopharyngeal Swab	Fever to be tested?	-	-	> 40	/
391	00051978	M	29	Nasopharyngeal Swab	Fever to be tested?	-	-	> 40	/
392	00063739	M	10	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/
393	00049026	M	36	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/
394	00057269	M	23	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/
395	00069359	F	34	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/
396	00010804	F	31	Nasopharyngeal Swab	Fever to be diagnosed	-	-	> 40	/
397	00057915	M	57	Nasopharyngeal Swab	Fever to be diagnosed, Enteritis?	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

398	00029877	M	69	Nasopharyngeal Swab	Lung Infection?	-	-	> 40	/
399	00081788	M	21	Nasopharyngeal Swab	Acute Tonsillitis	-	-	> 40	/
400	00063964	F	39	Nasopharyngeal Swab	URTI?	-	-	> 40	/
401	00582121	F	16	Nasopharyngeal Swab	URTI?	-	-	> 40	/
402	00036032	F	38	Nasopharyngeal Swab	Intrauterine adhesions following miscarriage, HIV	-	-	> 40	/
403	00050514	M	31	Nasopharyngeal Swab	Lung Infection, Fever	-	-	> 40	/
404	00059024	M	80	Nasopharyngeal Swab	Asthma to be tested?	-	-	> 40	/
405	00009268	M	73	Nasopharyngeal Swab	Community Acquired Pneumonia	-	-	> 40	/
406	00039601	F	47	Nasopharyngeal Swab	Chest tightness to be tested	-	-	> 40	/
407	00049014	M	61	Nasopharyngeal Swab	HIV, Liver Space Occupying: Primary Liver Cancer	-	-	> 40	/
408	00068238	F	54	Nasopharyngeal Swab	Lung Infection, Fever	-	-	> 40	/
409	00045269	M	59	Nasopharyngeal Swab	Back pain to be tested	-	-	> 40	/
410	00095475	M	88	Nasopharyngeal Swab	COPD with acute exacerbation	-	-	> 40	/
411	00657137	F	53	Nasopharyngeal Swab	Bronchial Asthma	-	-	> 40	/
412	00036708	M	64	Nasopharyngeal Swab	HBV (HBsAg, HBeAb, and HBcAb positive), Lung Infection	-	-	> 40	/
413	00069285	M	82	Nasopharyngeal Swab	Community Acquired Pneumonia, non-severe	-	-	> 40	/
414	00035612	M	65	Nasopharyngeal Swab	Foreign body in bronchus? Aspiration Pneumonia	-	-	> 40	/
415	00086439	M	83	Nasopharyngeal Swab	Community Acquired Pneumonia, non-severe	-	-	> 40	/
416	00044736	M	46	Nasopharyngeal Swab	Chronic HBV, fever	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

417	00026923	M	51	Nasopharyngeal Swab	Cough to be tested	-	-	> 40	/
418	00074580	F	77	Nasopharyngeal Swab	Renal Pelvis Tumor, Acute Gastroenteritis	-	-	> 40	/
419	00083392	F	30	Nasopharyngeal Swab	Ectopic Pregnancy	-	-	> 40	/
420	00032060	F	62	Nasopharyngeal Swab	HIV, Cervical Cancer	-	-	> 40	/

421	00017168	M	72	Nasopharyngeal Swab	Community Acquired Pneumonia	-	-	> 40	/
422	00047905	F	82	Nasopharyngeal Swab	COPD with acute exacerbation	-	-	> 40	/
423	00035241	F	88	Nasopharyngeal Swab	Fever to be tested? Severe Malnutrition	-	-	> 40	/
424	00059469	F	29	Nasopharyngeal Swab	Cough to be tested	-	-	> 40	/
425	00047932	F	34	Nasopharyngeal Swab	Community Acquired Pneumonia, non-severe	-	-	> 40	/
426	00056163	F	59	Nasopharyngeal Swab	Diabetes	-	-	> 40	/
427	00074525	F	36	Nasopharyngeal Swab	HBV (HBsAg, HBeAb, and HBcAb positive), Early Intrauterine Pregnancy	-	-	> 40	/
428	00042364	F	55	Nasopharyngeal Swab	Community Acquired Pneumonia, non-severe	-	-	> 40	/
429	00035253	M	84	Nasopharyngeal Swab	Lung Infection	-	-	> 40	/
430	00038920	M	63	Nasopharyngeal Swab	Pleural Effusion	-	-	> 40	/
431	00065262	F	50	Nasopharyngeal Swab	Acute Bronchitis? Community Acquired Pneumonia?	-	-	> 40	/
432	00088436	M	77	Nasopharyngeal Swab	Asthma to be tested	-	-	> 40	/
433	00042278	M	58	Nasopharyngeal Swab	COPD with acute exacerbation	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thüringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

434	00084568	F	92	Nasopharyngeal Swab	Skin ulcer of both legs to be tested	+	-	> 40	/
435	00025729	M	93	Nasopharyngeal Swab	COPD	-	-	> 40	/
436	00025884	M	83	Nasopharyngeal Swab	Lung Infection	-	-	> 40	/
437	00045779	F	41	Nasopharyngeal Swab	Missed Miscarriage?	-	-	> 40	/
438	00035568	M	91	Nasopharyngeal Swab	Fever to be tested, lung infection?	-	-	> 40	/
439	00073229	F	51	Nasopharyngeal Swab	Intrauterine Device Incarceration	-	-	> 40	/
440	00026516	F	29	Nasopharyngeal Swab	Expected Parturition	-	-	> 40	/
441	00083590	F	34	Nasopharyngeal Swab	Expected Parturition	-	-	> 40	/
442	00056577	M	50	Nasopharyngeal Swab	HBV (HBsAg, HBeAb, and HBeAb positive), Fever	-	-	> 40	/
443	00046085	M	64	Nasopharyngeal Swab	AIDS	-	-	> 40	/
444	00658638	M	27	Nasopharyngeal Swab	AIDS	-	-	> 40	/
445	00065360	M	33	Nasopharyngeal Swab	AIDS	-	-	> 40	/
446	00072834	M	54	Nasopharyngeal Swab	AIDS	-	-	> 40	/
447	00047965	M	78	Nasopharyngeal Swab	AIDS	-	-	> 40	/
448	00585637	M	30	Nasopharyngeal Swab	AIDS	-	-	> 40	/
449	00972618	M	48	Nasopharyngeal Swab	AIDS	-	-	> 40	/
450	00049724	F	42	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
451	00045436	M	42	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
452	00085957	M	34	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
453	00016659	F	44	Nasopharyngeal Swab	AIDS	-	-	> 40	/
454	00028359	M	26	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
455	00045398	M	52	Nasopharyngeal Swab	AIDS	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thüringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

456	00036615	F	29	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
457	00063164	M	27	Nasopharyngeal Swab	AIDS	-	-	> 40	/
458	00084627	M	34	Nasopharyngeal Swab	AIDS	-	-	> 40	/
459	QY005618	M	44	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
460	00035639	M	65	Nasopharyngeal Swab	AIDS	-	-	> 40	/
461	00654365	M	67	Nasopharyngeal Swab	Tuberculosis, HBV Carrier	-	-	> 40	/
462	00056918	M	30	Nasopharyngeal Swab	Secondary Tuberculosis (re-treatment, single drug resistance) smear- and culture-positive, Chronic HBV	-	-	> 40	/
463	00048454	M	40	Nasopharyngeal Swab	Tuberculosis, Lung Infection,	-	-	> 40	/
464	00069146	M	70	Nasopharyngeal Swab	Tuberculous Pleurisy	-	-	> 40	/
465	00062952	F	31	Nasopharyngeal Swab	Tuberculosis, Leukopenia, HBV	-	-	> 40	/
466	00075284	M	34	Nasopharyngeal Swab	Secondary Tuberculosis	-	-	> 40	/
467	00054267	F	58	Nasopharyngeal Swab	AIDS,	-	-	> 40	/
468	00048591	M	28	Nasopharyngeal Swab	HIV	-	-	> 40	/
469	00077425	M	64	Nasopharyngeal Swab	Tuberculosis recovered, off drugs and re-checked	-	-	> 40	/
470	00076157	M	33	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
471	00072258	M	57	Nasopharyngeal Swab	Tuberculous Meningitis	-	-	> 40	/
472	00081291	M	61	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
473	00051446	M	65	Nasopharyngeal Swab	Tuberculosis (drug resistance)	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

474	00063527	M	43	Nasopharyngeal Swab	AIDS	-	-	> 40	/
475	00016662	M	45	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
476	00038469	M	63	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
477	00082663	F	49	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
478	00017591	M	28	Nasopharyngeal Swab	Tuberculous Meningitis	-	-	> 40	/
479	00034652	F	22	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
480	00095655	F	55	Nasopharyngeal Swab	NTM	-	-	> 40	/
481	00063794	M	53	Nasopharyngeal Swab	AIDS	-	-	> 40	/
482	00075368	M	43	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
483	2003101740	F	72	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
484	2003129839	M	66	Nasopharyngeal Swab	Tuberculosis	-	-	> 40	/
485	2003173915	F	13	Nasopharyngeal Swab	Cervical Lymph Node Tuberculosis	-	-	> 40	/
486	2003122368	M	40	Nasopharyngeal Swab	Lung Shadow	-	-	> 40	/
487	2003195642	M	64	Nasopharyngeal Swab	Tuberculosis, Bronchial Tuberculosis	-	-	> 40	/
488	2003174628	F	65	Nasopharyngeal Swab	Chest Pain	-	-	> 40	/
489	2003060812	M	29	Nasopharyngeal Swab	AIDS	-	-	> 40	/
490	2003061512	F	14	Nasopharyngeal Swab	Bronchial Tuberculosis, Leukopenia	-	-	> 40	/
491	2003061841	M	58	Nasopharyngeal Swab	Chronic HBV	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

492	2003073731	M	51	Nasopharyngeal Swab	Secondary Tuberculosis (re-treatment, multiple drug resistance) smear-negative and culture-positive, Malnutrition, Type II Diabetes	-	-	> 40	/
493	2003097466	M	28	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
494	2003099731	M	33	Nasopharyngeal Swab	Sore throat to be tested	-	-	> 40	/
495	2003114417	M	37	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
496	2003114831	F	26	Nasopharyngeal Swab	To be tested	-	-	> 40	/
497	2003101628	M	42	Nasopharyngeal Swab	Fever, diarrhea, abdominal distention to be tested	-	-	> 40	/
498	2003102627	M	22	Nasopharyngeal Swab	Lung lesions to be tested	-	-	> 40	/
499	200311573	M	22	Nasopharyngeal Swab	Fever, chest pain	-	-	> 40	/
500	2003116517	F	23	Nasopharyngeal Swab	Fever to be tested (URTI? Salmonellosis?)	-	-	> 40	/
501	2003128163	M	37	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
502	2003128736	M	43	Nasopharyngeal Swab	Chest pain to be tested	-	-	> 40	/
503	2003129368	F	28	Nasopharyngeal Swab	Pyogenic Tonsillitis	-	-	> 40	/
504	2003120471	M	49	Nasopharyngeal Swab	Pharyngeal Bronchitis	-	-	> 40	/
505	2003133958	M	18	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
506	2003145474	M	22	Nasopharyngeal Swab	Fever and emesis to be tested	-	-	> 40	/
507	2003161259	M	24	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
508	2003173956	M	48	Nasopharyngeal Swab	Chest tightness to be tested	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

509	2003174854	M	38	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
510	2003174955	F	48	Nasopharyngeal Swab	Fever: Pharynx Tonsillitis	-	-	> 40	/
511	2003191633	M	85	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
512	202007340	M	18	Nasopharyngeal Swab	Cough to be tested	-	-	> 40	/
513	202008570	M	36	Nasopharyngeal Swab	Pyogenic Tonsillitis	-	-	> 40	/
514	202006479	F	23	Nasopharyngeal Swab	Fever, Sore Throat	-	-	> 40	/
515	202008578	M	30	Nasopharyngeal Swab	Fever to be tested	-	-	> 40	/
516	202008759	F	47	Nasopharyngeal Swab	Chronic HBV with δ factor	-	-	> 40	/
517	202008236	M	66	Nasopharyngeal Swab	Thoracic Abscess	-	-	> 40	/
518	201104001	M	57	Nasopharyngeal Swab	/	-	-	> 40	/
519	201104002	M	53	Nasopharyngeal Swab	/	-	-	> 40	/
520	201104003	M	43	Nasopharyngeal Swab	/	-	-	> 40	/
521	201104004	F	44	Nasopharyngeal Swab	/	-	-	> 40	/
522	201104005	F	39	Nasopharyngeal Swab	/	-	-	> 40	/
523	201104006	F	26	Nasopharyngeal Swab	/	-	-	> 40	/
524	201104007	M	56	Nasopharyngeal Swab	/	-	-	> 40	/
525	201104008	F	30	Nasopharyngeal Swab	/	-	-	> 40	/
526	201104009	F	45	Nasopharyngeal Swab	/	-	-	> 40	/
527	201104010	F	40	Nasopharyngeal Swab	/	-	-	> 40	/
528	201104011	F	35	Nasopharyngeal Swab	/	-	-	> 40	/
529	201104012	F	29	Nasopharyngeal Swab	/	-	-	> 40	/
530	201104013	F	22	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

531	201104014	F	43	Nasopharyngeal Swab	/	-	-	> 40	/
532	201104015	F	45	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
533	201104016	F	45	Nasopharyngeal Swab	/	-	-	> 40	/
534	201104017	F	43	Nasopharyngeal Swab	/	-	-	> 40	/
535	201104018	M	43	Nasopharyngeal Swab	/	-	-	> 40	/
536	201104019	M	41	Nasopharyngeal Swab	/	-	-	> 40	/
537	201104020	M	33	Nasopharyngeal Swab	/	-	-	> 40	/

538	201104021	M	32	Nasopharyngeal Swab	/	-	-	> 40	/
539	201104022	F	61	Nasopharyngeal Swab	/	-	-	> 40	/
540	201104023	F	44	Nasopharyngeal Swab	/	-	-	> 40	/
541	201104024	M	37	Nasopharyngeal Swab	/	-	-	> 40	/
542	201104025	M	28	Nasopharyngeal Swab	/	-	-	> 40	/
543	201104026	M	37	Nasopharyngeal Swab	Cold	-	-	> 40	/
544	201104027	M	29	Nasopharyngeal Swab	/	-	-	> 40	/
545	201104028	M	29	Nasopharyngeal Swab	/	-	-	> 40	/
546	201104029	M	29	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
547	201104030	M	36	Nasopharyngeal Swab	/	-	-	> 40	/
548	201104031	M	26	Nasopharyngeal Swab	/	-	-	> 40	/
549	201104032	M	26	Nasopharyngeal Swab	/	-	-	> 40	/
550	201104033	M	33	Nasopharyngeal Swab	/	-	-	> 40	/
551	201104034	M	24	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

552	201104035	M	25	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
553	201104036	M	24	Nasopharyngeal Swab	/	-	-	> 40	/
554	201104037	M	24	Nasopharyngeal Swab	/	-	-	> 40	/
555	201104038	M	27	Nasopharyngeal Swab	/	-	-	> 40	/
556	201104039	M	24	Nasopharyngeal Swab	/	-	-	> 40	/
557	201104040	M	23	Nasopharyngeal Swab	/	-	-	> 40	/
558	201104041	M	32	Nasopharyngeal Swab	/	-	-	> 40	/
559	201104042	M	25	Nasopharyngeal Swab	/	-	-	> 40	/
560	201104043	M	26	Nasopharyngeal Swab	/	-	-	> 40	/
561	201104044	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
562	201104045	F	54	Nasopharyngeal Swab	/	-	-	> 40	/
563	201104046	M	31	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
564	201104047	F	32	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
565	201104048	M	49	Nasopharyngeal Swab	/	-	-	> 40	/
566	201104049	M	52	Nasopharyngeal Swab	/	-	-	> 40	/
567	201104050	M	42	Nasopharyngeal Swab	/	-	-	> 40	/
568	201104051	M	48	Nasopharyngeal Swab	/	-	-	> 40	/
569	201104052	M	39	Nasopharyngeal Swab	/	-	-	> 40	/
570	201104053	M	39	Nasopharyngeal Swab	/	-	-	> 40	/
571	201104054	F	42	Nasopharyngeal Swab	/	-	-	> 40	/
572	201104055	M	30	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

573	201104056	M	47	Nasopharyngeal Swab	/	-	-	> 40	/
574	201104057	M	38	Nasopharyngeal Swab	/	-	-	> 40	/
575	201104058	M	49	Nasopharyngeal Swab	/	-	-	> 40	/
576	201104059	M	27	Nasopharyngeal Swab	/	-	-	> 40	/
577	201104060	M	45	Nasopharyngeal Swab	/	-	-	> 40	/
578	201104061	M	40	Nasopharyngeal Swab	/	-	-	> 40	/
579	201104062	M	38	Nasopharyngeal Swab	/	-	-	> 40	/
580	201104063	M	38	Nasopharyngeal Swab	/	-	-	> 40	/
581	201104064	F	40	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
582	201104065	F	31	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
583	201104066	F	27	Nasopharyngeal Swab	/	-	-	> 40	/
584	201104067	F	29	Nasopharyngeal Swab	/	-	-	> 40	/
585	201104068	F	39	Nasopharyngeal Swab	/	-	-	> 40	/
586	201104069	F	32	Nasopharyngeal Swab	/	-	-	> 40	/
587	201104070	F	39	Nasopharyngeal Swab	/	-	-	> 40	/
588	201104071	F	42	Nasopharyngeal Swab	/	-	-	> 40	/
589	201104072	F	36	Nasopharyngeal Swab	/	-	-	> 40	/
590	201104073	F	31	Nasopharyngeal Swab	Cold	-	-	> 40	/
591	201104074	F	30	Nasopharyngeal Swab	/	-	-	> 40	/
592	201104075	F	35	Nasopharyngeal Swab	/	-	-	> 40	/
593	201104076	F	27	Nasopharyngeal Swab	/	-	-	> 40	/
594	201104077	F	28	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

595	201104078	M	24	Nasopharyngeal Swab	/	-	-	> 40	/
596	201104079	M	28	Nasopharyngeal Swab	/	-	-	> 40	/
597	201104080	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
598	201104081	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
599	201104082	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
600	201104083	F	27	Nasopharyngeal Swab	/	-	-	> 40	/
601	201104084	F	37	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
602	201104085	F	26	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
603	201104086	M	29	Nasopharyngeal Swab	/	-	-	> 40	/
604	201104087	M	36	Nasopharyngeal Swab	/	-	-	> 40	/
605	201104088	F	28	Nasopharyngeal Swab	/	-	-	> 40	/
606	201104089	F	26	Nasopharyngeal Swab	/	-	-	> 40	/
607	201104090	M	23	Nasopharyngeal Swab	/	-	-	> 40	/
608	201104091	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
609	201104092	F	55	Nasopharyngeal Swab	/	-	-	> 40	/
610	201104093	M	66	Nasopharyngeal Swab	/	-	-	> 40	/
611	201104094	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
612	201104095	F	28	Nasopharyngeal Swab	/	-	-	> 40	/
613	201104096	F	27	Nasopharyngeal Swab	/	-	-	> 40	/
614	201104097	F	25	Nasopharyngeal Swab	/	-	-	> 40	/
615	201104098	F	25	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

616	201104099	M	46	Nasopharyngeal Swab	/	-	-	> 40	/
617	201104100	M	54	Nasopharyngeal Swab	/	-	-	> 40	/
618	201104101	M	25	Nasopharyngeal Swab	/	-	-	> 40	/
619	201104102	F	25	Nasopharyngeal Swab	/	-	-	> 40	/
620	201104103	F	22	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
621	201104104	M	30	Nasopharyngeal Swab	/	-	-	> 40	/
622	201104105	F	34	Nasopharyngeal Swab	/	-	-	> 40	/
623	201104106	F	29	Nasopharyngeal Swab	/	-	-	> 40	/
624	201104107	F	27	Nasopharyngeal Swab	/	-	-	> 40	/

625	201104108	F	29	Nasopharyngeal Swab	/	-	-	> 40	/
626	201104109	F	32	Nasopharyngeal Swab	/	-	-	> 40	/
627	201104110	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
628	201104111	F	26	Nasopharyngeal Swab	/	-	-	> 40	/
629	201104112	M	25	Nasopharyngeal Swab	/	-	-	> 40	/
630	201104113	F	29	Nasopharyngeal Swab	/	-	-	> 40	/
631	201104114	F	26	Nasopharyngeal Swab	/	-	-	> 40	/
632	201104115	M	28	Nasopharyngeal Swab	/	-	-	> 40	/
633	201104116	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
634	201104117	M	27	Nasopharyngeal Swab	/	-	-	> 40	/
635	201104118	F	30	Nasopharyngeal Swab	/	-	-	> 40	/
636	201104119	M	23	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

637	201104120	M	25	Nasopharyngeal Swab	/	-	-	> 40	/
638	201104121	M	30	Nasopharyngeal Swab	/	-	-	> 40	/
639	201104122	M	25	Nasopharyngeal Swab	/	-	-	> 40	/
640	201104123	M	27	Nasopharyngeal Swab	/	-	-	> 40	/
641	201104124	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
642	201104125	F	25	Nasopharyngeal Swab	/	-	-	> 40	/
643	201104126	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
644	201104127	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
645	201104128	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
646	201104129	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
647	201104130	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
648	201104131	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
649	201104132	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
650	201104133	M	24	Nasopharyngeal Swab	/	-	-	> 40	/
651	201104134	M	24	Nasopharyngeal Swab	/	-	-	> 40	/
652	201104135	M	31	Nasopharyngeal Swab	/	-	-	> 40	/
653	201104136	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
654	201104137	M	35	Nasopharyngeal Swab	Pharyngitis	-	-	> 40	/
655	201104138	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
656	201104139	F	27	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
657	201104140	M	22	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

658	201104141	M	27	Nasopharyngeal Swab	/	-	-	> 40	/
659	201104142	F	33	Nasopharyngeal Swab	/	-	-	> 40	/
660	201104143	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
661	201104144	M	30	Nasopharyngeal Swab	/	-	-	> 40	/
662	201104145	M	32	Nasopharyngeal Swab	/	-	-	> 40	/
663	201104146	M	26	Nasopharyngeal Swab	/	-	-	> 40	/
664	201104147	F	26	Nasopharyngeal Swab	/	-	-	> 40	/
665	201104148	F	25	Nasopharyngeal Swab	/	-	-	> 40	/
666	201104149	F	25	Nasopharyngeal Swab	/	-	-	> 40	/
667	201104150	M	31	Nasopharyngeal Swab	/	-	-	> 40	/
668	201104151	F	33	Nasopharyngeal Swab	/	-	-	> 40	/
669	201104152	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
670	201104153	F	25	Nasopharyngeal Swab	/	-	-	> 40	/
671	201104154	F	27	Nasopharyngeal Swab	/	-	-	> 40	/
672	201104155	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
673	201104156	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
674	201104157	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
675	201104158	M	28	Nasopharyngeal Swab	/	-	-	> 40	/
676	201104159	M	37	Nasopharyngeal Swab	/	-	-	> 40	/
677	201104160	F	36	Nasopharyngeal Swab	/	-	-	> 40	/
678	201104161	F	31	Nasopharyngeal Swab	Pharyngitis	-	-	> 40	/
679	201104162	F	27	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

680	201104163	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
681	201104164	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
682	201104165	F	24	Nasopharyngeal Swab	/	-	-	> 40	/

683	201104166	F	31	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
684	201104167	F	43	Nasopharyngeal Swab	/	-	-	> 40	/
685	201104168	F	33	Nasopharyngeal Swab	/	-	-	> 40	/
686	201104169	F	26	Nasopharyngeal Swab	/	-	-	> 40	/
687	201104170	F	25	Nasopharyngeal Swab	/	-	-	> 40	/
688	201104171	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
689	201104172	F	26	Nasopharyngeal Swab	/	-	-	> 40	/
690	201104173	F	33	Nasopharyngeal Swab	/	-	-	> 40	/
691	201104174	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
692	201104175	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
693	201104176	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
694	201104177	F	32	Nasopharyngeal Swab	/	-	-	> 40	/
695	201104178	M	35	Nasopharyngeal Swab	/	-	-	> 40	/
696	201104179	M	30	Nasopharyngeal Swab	/	-	-	> 40	/
697	201104180	M	39	Nasopharyngeal Swab	/	-	-	> 40	/
698	201104181	M	27	Nasopharyngeal Swab	/	-	-	> 40	/
699	201104182	M	54	Nasopharyngeal Swab	/	-	-	> 40	/
700	201104183	M	49	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thüringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

701	201104184	M	39	Nasopharyngeal Swab	/	-	-	> 40	/
702	201104185	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
703	201104186	F	38	Nasopharyngeal Swab	Rhinitis	-	-	> 40	/
704	201104187	F	36	Nasopharyngeal Swab	/	-	-	> 40	/
705	201104188	F	27	Nasopharyngeal Swab	/	-	-	> 40	/
706	201104189	F	28	Nasopharyngeal Swab	/	-	-	> 40	/
707	201104190	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
708	201104191	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
709	201104192	F	22	Nasopharyngeal Swab	/	-	-	> 40	/
710	201104193	F	24	Nasopharyngeal Swab	/	-	-	> 40	/
711	201104194	F	23	Nasopharyngeal Swab	/	-	-	> 40	/

712	201104195	F	23	Nasopharyngeal Swab	/	-	-	> 40	/
713	201104196	F	28	Nasopharyngeal Swab	/	-	-	> 40	/
714	201104197	F	38	Nasopharyngeal Swab	/	-	-	> 40	/
715	201104198	F	38	Nasopharyngeal Swab	/	-	-	> 40	/
716	201104199	F	32	Nasopharyngeal Swab	/	-	-	> 40	/
717	201104200	F	42	Nasopharyngeal Swab	/	-	-	> 40	/
718	201104201	M	67	Nasopharyngeal Swab	/	-	-	> 40	/
719	201104202	F	41	Nasopharyngeal Swab	/	-	-	> 40	/
720	201104203	F	38	Nasopharyngeal Swab	/	-	-	> 40	/
721	201104204	F	45	Nasopharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

722	201104205	F	39	Nasopharyngeal Swab	/	-	-	> 40	/
723	201104206	F	31	Nasopharyngeal Swab	/	-	-	> 40	/
724	201104207	F	36	Nasopharyngeal Swab	/	-	-	> 40	/
725	201104208	F	37	Nasopharyngeal Swab	/	-	-	> 40	/
726	201104209	F	39	Nasopharyngeal Swab	/	-	-	> 40	/
727	201104210	F	54	Nasopharyngeal Swab	/	-	-	> 40	/
728	201202035	F	79	Nasopharyngeal Swab	/	-	-	> 35	/
729	201202036	M	68	Nasopharyngeal Swab	/	-	-	> 35	/
730	201202037	M	55	Nasopharyngeal Swab	/	-	-	> 35	/
731	201202038	F	39	Nasopharyngeal Swab	/	-	-	> 35	/
732	201202039	F	52	Nasopharyngeal Swab	/	-	-	> 35	/
733	201202040	F	35	Nasopharyngeal Swab	/	-	-	> 35	/
734	201202041	F	57	Nasopharyngeal Swab	/	-	-	> 35	/
735	201202042	M	62	Nasopharyngeal Swab	/	-	-	> 35	/
736	201202043	F	51	Nasopharyngeal Swab	/	-	-	> 35	/
737	201202044	M	51	Nasopharyngeal Swab	/	-	-	> 35	/
738	201202045	F	84	Nasopharyngeal Swab	/	-	-	> 35	/
739	201202046	F	74	Nasopharyngeal Swab	/	-	-	> 35	/
740	201202047	F	47	Nasopharyngeal Swab	/	-	-	> 35	/
741	201202048	M	59	Nasopharyngeal Swab	/	-	-	> 35	/
742	201202049	F	85	Nasopharyngeal Swab	/	-	-	> 35	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

743	201202050	M	60	Nasopharyngeal Swab	/	-	-	> 35	/
744	201202051	F	55	Nasopharyngeal Swab	/	-	-	> 35	/
745	201202052	M	18	Nasopharyngeal Swab	/	-	-	> 35	/
746	201202053	F	25	Nasopharyngeal Swab	/	-	-	> 35	/
747	201202054	F	40	Nasopharyngeal Swab	/	-	-	> 35	/
748	201202055	F	82	Nasopharyngeal Swab	/	-	-	> 35	/
749	201202056	F	89	Nasopharyngeal Swab	/	-	-	> 35	/
750	201202057	M	56	Nasopharyngeal Swab	/	-	-	> 35	/
751	201202058	M	62	Nasopharyngeal Swab	/	-	-	> 35	/
752	201202059	F	46	Nasopharyngeal Swab	/	-	-	> 35	/
753	201202060	F	8	Nasopharyngeal Swab	/	-	-	> 35	/
754	201202061	F	42	Nasopharyngeal Swab	/	-	-	> 35	/
755	201202062	F	25	Nasopharyngeal Swab	/	-	-	> 35	/
756	201202063	M	33	Nasopharyngeal Swab	/	-	-	> 35	/
757	201202064	F	15	Nasopharyngeal Swab	/	-	-	> 35	/
758	201202065	M	25	Nasopharyngeal Swab	/	-	-	> 35	/
759	201202066	F	43	Nasopharyngeal Swab	/	-	-	> 35	/
760	201202067	M	79	Nasopharyngeal Swab	/	-	-	> 35	/
761	201202068	F	47	Nasopharyngeal Swab	/	-	-	> 35	/
762	201202069	M	51	Nasopharyngeal Swab	/	-	-	> 35	/
763	201202070	F	76	Nasopharyngeal Swab	/	-	-	> 35	/
764	201202071	M	48	Nasopharyngeal Swab	/	-	-	> 35	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

765	201202072	F	82	Nasopharyngeal Swab	/	-	-	> 35	/
766	201202073	F	50	Nasopharyngeal Swab	/	-	-	> 35	/
767	201202074	F	66	Nasopharyngeal Swab	/	-	-	> 35	/
768	201202075	F	84	Nasopharyngeal Swab	/	-	-	> 35	/
769	201202076	M	22	Nasopharyngeal Swab	/	-	-	> 35	/

770	201202077	F	68	Nasopharyngeal Swab	/	-	-	> 35	/
771	201202078	M	24	Nasopharyngeal Swab	/	-	-	> 35	/
772	201202079	F	70	Nasopharyngeal Swab	/	-	-	> 35	/
773	201202080	M	29	Nasopharyngeal Swab	/	-	-	> 35	/
774	201202081	M	76	Nasopharyngeal Swab	/	-	-	> 35	/
775	201202082	F	46	Nasopharyngeal Swab	/	-	-	> 35	/
776	201202083	F	62	Nasopharyngeal Swab	/	-	-	> 35	/
777	201202084	F	61	Nasopharyngeal Swab	/	-	-	> 35	/
778	210112001	F	23	Nasal Swab	/	-	-	> 40	/
779	210112002	F	22	Nasal Swab	/	-	-	> 40	/
780	210112003	M	28	Nasal Swab	/	-	-	> 40	/
781	210112004	M	30	Nasal Swab	/	-	-	> 40	/
782	210112005	F	21	Nasal Swab	/	-	-	> 40	/
783	210112006	F	21	Nasal Swab	/	-	-	> 40	/
784	210112007	M	34	Nasal Swab	/	-	-	> 40	/
785	210112008	F	37	Nasal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

786	210112009	M	31	Nasal Swab	/	-	-	> 40	/
787	210112010	M	31	Nasal Swab	/	-	-	> 40	/
788	210112011	F	27	Nasal Swab	/	-	-	> 40	/
789	210112012	F	35	Nasal Swab	/	-	-	> 40	/
790	210112013	M	45	Nasal Swab	/	-	-	> 40	/
791	210112014	F	30	Nasal Swab	/	-	-	> 40	/
792	210112015	M	49	Nasal Swab	/	-	-	> 40	/
793	210112016	F	38	Nasal Swab	/	-	-	> 40	/
794	210112017	F	43	Nasal Swab	/	+	-	> 40	/
795	210112018	F	36	Nasal Swab	/	-	-	> 40	/
796	210112019	F	19	Nasal Swab	/	-	-	> 40	/
797	210112020	M	46	Nasal Swab	/	-	-	> 40	/
798	210112021	F	40	Nasal Swab	/	-	-	> 40	/

799	210112022	M	24	Nasal Swab	/	-	-	> 40	/
800	210112023	F	44	Nasal Swab	/	-	-	> 40	/
801	210112024	F	38	Nasal Swab	/	-	-	> 40	/
802	210112025	M	31	Nasal Swab	/	-	-	> 40	/
803	210112026	F	26	Nasal Swab	/	-	-	> 40	/
804	210112027	F	24	Nasal Swab	/	-	-	> 40	/
805	210112028	F	31	Nasal Swab	/	-	-	> 40	/
806	210112029	M	29	Nasal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

807	210112030	M	36	Nasal Swab	/	-	-	> 40	/
808	210112031	F	40	Nasal Swab	/	-	-	> 40	/
809	210112032	F	45	Nasal Swab	/	-	-	> 40	/
810	210112033	M	31	Nasal Swab	/	-	-	> 40	/
811	210112034	F	35	Nasal Swab	/	-	-	> 40	/
812	210112035	F	29	Nasal Swab	/	-	-	> 40	/
813	210112036	F	49	Nasal Swab	/	-	-	> 40	/
814	210112037	M	41	Nasal Swab	/	-	-	> 40	/
815	210112038	F	25	Nasal Swab	/	-	-	> 40	/
816	210112039	F	33	Nasal Swab	/	-	-	> 40	/
817	210112040	F	30	Nasal Swab	/	-	-	> 40	/
818	210112041	F	23	Nasal Swab	/	-	-	> 40	/
819	210112042	M	34	Nasal Swab	/	-	-	> 40	/
820	210112043	M	31	Nasal Swab	/	-	-	> 40	/
821	210112044	F	27	Nasal Swab	/	-	-	> 40	/
822	210112045	F	32	Nasal Swab	/	-	-	> 40	/
823	210112046	F	46	Nasal Swab	/	-	-	> 40	/
824	210112047	F	47	Nasal Swab	/	-	-	> 40	/
825	210112048	F	25	Nasal Swab	/	-	-	> 40	/
826	210112049	M	33	Nasal Swab	/	-	-	> 40	/
827	210112050	F	38	Nasal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

828	210112051	M	30	Nasal Swab	/	-	-	> 40	/
829	210112052	M	38	Nasal Swab	/	-	-	> 40	/
830	210112053	F	31	Nasal Swab	/	-	-	> 40	/
831	210112054	M	38	Nasal Swab	/	-	-	> 40	/
832	210112055	M	31	Nasal Swab	/	-	-	> 40	/
833	210112056	M	24	Nasal Swab	/	-	-	> 40	/
834	210112057	F	35	Nasal Swab	/	-	-	> 40	/
835	210112058	F	55	Nasal Swab	/	-	-	> 40	/
836	210112059	M	39	Nasal Swab	/	-	-	> 40	/
837	210112060	F	42	Nasal Swab	/	-	-	> 40	/
838	210112061	M	37	Nasal Swab	/	-	-	> 40	/
839	210112062	M	31	Nasal Swab	/	-	-	> 40	/
840	210112063	M	45	Nasal Swab	/	-	-	> 40	/
841	210112064	F	30	Nasal Swab	/	-	-	> 40	/
842	210112065	F	45	Nasal Swab	/	-	-	> 40	/
843	210112066	M	47	Nasal Swab	/	-	-	> 40	/
844	210112067	F	37	Nasal Swab	/	-	-	> 40	/
845	210112068	F	25	Nasal Swab	/	-	-	> 40	/
846	210112069	F	51	Nasal Swab	/	-	-	> 40	/
847	210112070	M	42	Nasal Swab	/	-	-	> 40	/
848	210112071	F	39	Nasal Swab	/	-	-	> 40	/
849	210112072	F	25	Nasal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

850	210112073	M	20	Nasal Swab	/	-	-	> 40	/
851	210112074	F	23	Nasal Swab	/	-	-	> 40	/
852	210112075		29	Nasal Swab	/	-	-	> 40	/
853	210112076	M	23	Nasal Swab	/	-	-	> 40	/
854	210112077	F	46	Nasal Swab	/	-	-	> 40	/
855	210112078	M	27	Nasal Swab	/	-	-	> 40	/
856	210112079	M	30	Nasal Swab	/	-	-	> 40	/

857	210112080	F	24	Nasal Swab	/	-	-	> 40	/
858	210112081	F	49	Nasal Swab	/	-	-	> 40	/
859	210112082	F	26	Nasal Swab	/	-	-	> 40	/
860	210112083	M	29	Nasal Swab	/	-	-	> 40	/
861	210112084	M	33	Nasal Swab	/	-	-	> 40	/
862	210112085	F	24	Nasal Swab	/	-	-	> 40	/
863	210112086	F	24	Nasal Swab	/	-	-	> 40	/
864	210112087	M	31	Nasal Swab	/	-	-	> 40	/
865	210112088	M	28	Nasal Swab	/	-	-	> 40	/
866	210112089	F	26	Nasal Swab	/	-	-	> 40	/
867	210112090	F	38	Nasal Swab	/	-	-	> 40	/
868	210112091	F	40	Nasal Swab	/	-	-	> 40	/
869	210112092	F	38	Nasal Swab	/	-	-	> 40	/
870	210112093	F	35	Nasal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

871	210112094	M	23	Nasal Swab	/	-	-	> 40	/
872	210112095	M	40	Nasal Swab	/	-	-	> 40	/
873	210112096	M	45	Nasal Swab	/	-	-	> 40	/
874	210112097	M	25	Nasal Swab	/	-	-	> 40	/
875	210112098	M	51	Nasal Swab	/	-	-	> 40	/
876	210112099	F	49	Nasal Swab	/	-	-	> 40	/
877	210112100	F	50	Nasal Swab	/	-	-	> 40	/
878	210112101	F	46	Nasal Swab	/	-	-	> 40	/
879	210112102	M	30	Nasal Swab	/	-	-	> 40	/
880	210112103	M	31	Nasal Swab	/	-	-	> 40	/
881	210112104	M	38	Nasal Swab	/	-	-	> 40	/
882	210112105	F	40	Nasal Swab	/	-	-	> 40	/
883	210112106	M	22	Nasal Swab	/	-	-	> 40	/
884	210112107	M	31	Nasal Swab	/	-	-	> 40	/
885	210112108	F	31	Nasal Swab	/	-	-	> 40	/
886	210112109	F	46	Nasal Swab	/	-	-	> 40	/
887	210112110	F	29	Nasal Swab	/	-	-	> 40	/
888	210112111	M	45	Nasal Swab	/	-	-	> 40	/
889	210112112	M	28	Nasal Swab	/	-	-	> 40	/
890	210112113	F	27	Nasal Swab	/	-	-	> 40	/
891	210112114	M	52	Nasal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

892	210112115	M	27	Nasal Swab	/	-	-	> 40	/
893	210112116	M	24	Nasal Swab	/	-	-	> 40	/
894	210112117	F	51	Nasal Swab	/	-	-	> 40	/
895	210112118	M	25	Nasal Swab	/	-	-	> 40	/
896	210112119	F	29	Nasal Swab	/	-	-	> 40	/
897	210112120	F	26	Nasal Swab	/	-	-	> 40	/
898	210112121	M	33	Nasal Swab	/	-	-	> 40	/
899	210112122	M	37	Nasal Swab	/	-	-	> 40	/
900	210112123	F	35	Nasal Swab	/	-	-	> 40	/
901	210112124	F	32	Nasal Swab	/	-	-	> 40	/
902	210112125	F	47	Nasal Swab	/	-	-	> 40	/
903	210112126	M	46	Oropharyngeal Swab	/	-	-	> 40	/
904	210112127	F	42	Oropharyngeal Swab	/	-	-	> 40	/
905	210112128	M	26	Oropharyngeal Swab	/	-	-	> 40	/
906	210112129	F	31	Oropharyngeal Swab	/	-	-	> 40	/
907	210112130	F	55	Oropharyngeal Swab	/	-	-	> 40	/
908	210112131	M	27	Oropharyngeal Swab	/	-	-	> 40	/
909	210112132	F	43	Oropharyngeal Swab	/	-	-	> 40	/
910	210112133	M	31	Oropharyngeal Swab	/	-	-	> 40	/
911	210112134	F	30	Oropharyngeal Swab	/	-	-	> 40	/
912	210112135	F	48	Oropharyngeal Swab	/	-	-	> 40	/
913	210112136	F	24	Oropharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

914	210112137	F	38	Oropharyngeal Swab	/	-	-	> 40	/
915	210112138	F	45	Oropharyngeal Swab	/	-	-	> 40	/
916	210112139	M	32	Oropharyngeal Swab	/	-	-	> 40	/
917	210112140	M	39	Oropharyngeal Swab	/	-	-	> 40	/
918	210112141	M	27	Oropharyngeal Swab	/	-	-	> 40	/
919	210112142	F	26	Oropharyngeal Swab	/	-	-	> 40	/
920	210112143	F	52	Oropharyngeal Swab	/	-	-	> 40	/
921	210112144	F	47	Oropharyngeal Swab	/	-	-	> 40	/
922	210112145	F	24	Oropharyngeal Swab	/	-	-	> 40	/
923	210112146	F	27	Oropharyngeal Swab	/	-	-	> 40	/
924	210112147	F	22	Oropharyngeal Swab	/	-	-	> 40	/
925	210112148	M	51	Oropharyngeal Swab	/	-	-	> 40	/
926	210112149	M	23	Oropharyngeal Swab	/	-	-	> 40	/
927	210112150	M	35	Oropharyngeal Swab	/	-	-	> 40	/
928	210112151	M	47	Oropharyngeal Swab	/	-	-	> 40	/
929	210112152	F	24	Oropharyngeal Swab	/	-	-	> 40	/
930	210112153	M	25	Oropharyngeal Swab	/	-	-	> 40	/
931	210112154	M	31	Oropharyngeal Swab	/	-	-	> 40	/
932	210112155	F	27	Oropharyngeal Swab	/	-	-	> 40	/
933	210112156	F	26	Oropharyngeal Swab	/	-	-	> 40	/
934	210112157	F	24	Oropharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

935	210112158	F	44	Oropharyngeal Swab	/	-	-	> 40	/
936	210112159	F	39	Oropharyngeal Swab	/	-	-	> 40	/
937	210112160	M	38	Oropharyngeal Swab	/	-	-	> 40	/
938	210112161	F	31	Oropharyngeal Swab	/	-	-	> 40	/
939	210112162	F	45	Oropharyngeal Swab	/	-	-	> 40	/
940	210112163	M	28	Oropharyngeal Swab	/	-	-	> 40	/
941	210112164	M	21	Oropharyngeal Swab	/	-	-	> 40	/
942	210112165	M	46	Oropharyngeal Swab	/	-	-	> 40	/
943	210112166	F	24	Oropharyngeal Swab	/	-	-	> 40	/

944	210112167	F	25	Oropharyngeal Swab	/	-	-	> 40	/
945	210112168	F	39	Oropharyngeal Swab	/	-	-	> 40	/
946	210112169	M	27	Oropharyngeal Swab	/	-	-	> 40	/
947	210112170	F	29	Oropharyngeal Swab	/	-	-	> 40	/
948	210112171	F	28	Oropharyngeal Swab	/	-	-	> 40	/
949	210112172	M	25	Oropharyngeal Swab	/	-	-	> 40	/
950	210112173	F	38	Oropharyngeal Swab	/	-	-	> 40	/
951	210112174	M	31	Oropharyngeal Swab	/	-	-	> 40	/
952	210112175	M	35	Oropharyngeal Swab	/	-	-	> 40	/
953	210112176	M	31	Oropharyngeal Swab	/	-	-	> 40	/
954	210112177	M	36	Oropharyngeal Swab	/	-	-	> 40	/
955	210112178	F	33	Oropharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

956	210112179	F	46	Oropharyngeal Swab	/	-	-	> 40	/
957	210112180	F	33	Oropharyngeal Swab	/	-	-	> 40	/
958	210112181	F	31	Oropharyngeal Swab	/	-	-	> 40	/
959	210112182	F	42	Oropharyngeal Swab	/	-	-	> 40	/
960	210112183	F	42	Oropharyngeal Swab	/	-	-	> 40	/
961	210112184	F	37	Oropharyngeal Swab	/	-	-	> 40	/
962	210112185	F	39	Oropharyngeal Swab	/	-	-	> 40	/
963	210112186	M	44	Oropharyngeal Swab	/	-	-	> 40	/
964	210112187	M	33	Oropharyngeal Swab	/	-	-	> 40	/
965	210112188	F	30	Oropharyngeal Swab	/	-	-	> 40	/
966	210112189	M	35	Oropharyngeal Swab	/	-	-	> 40	/
967	210112190	M	50	Oropharyngeal Swab	/	-	-	> 40	/
968	210112191	F	46	Oropharyngeal Swab	/	-	-	> 40	/
969	210112192	M	16	Oropharyngeal Swab	/	-	-	> 40	/
970	210112193	F	30	Oropharyngeal Swab	/	-	-	> 40	/
971	210112194	F	39	Oropharyngeal Swab	/	-	-	> 40	/
972	210112195	F	40	Oropharyngeal Swab	/	-	-	> 40	/
973	210112196	F	20	Oropharyngeal Swab	/	-	-	> 40	/
974	210112197	M	28	Oropharyngeal Swab	/	-	-	> 40	/
975	210112198	M	22	Oropharyngeal Swab	/	-	-	> 40	/
976	210112199	M	36	Oropharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

977	210112200	M	48	Oropharyngeal Swab	/	-	-	> 40	/
978	210112201	F	38	Oropharyngeal Swab	/	-	-	> 40	/
979	210112202	M	46	Oropharyngeal Swab	/	-	-	> 40	/
980	210112203	M	26	Oropharyngeal Swab	/	-	-	> 40	/
981	210112204	F	38	Oropharyngeal Swab	/	-	-	> 40	/
982	210112205	M	36	Oropharyngeal Swab	/	-	-	> 40	/
983	210112206	M	33	Oropharyngeal Swab	/	-	-	> 40	/
984	210112207	M	31	Oropharyngeal Swab	/	-	-	> 40	/
985	210112208	F	33	Oropharyngeal Swab	/	-	-	> 40	/
986	210112209	M	32	Oropharyngeal Swab	/	-	-	> 40	/
987	210112210	M	20	Oropharyngeal Swab	/	-	-	> 40	/
988	210112211	M	20	Oropharyngeal Swab	/	-	-	> 40	/
989	210112212	M	20	Oropharyngeal Swab	/	-	-	> 40	/
990	210112213	M	26	Oropharyngeal Swab	/	-	-	> 40	/
991	210112214	F	35	Oropharyngeal Swab	/	-	-	> 40	/
992	210112215	F	24	Oropharyngeal Swab	/	-	-	> 40	/
993	210112216	M	37	Oropharyngeal Swab	/	-	-	> 40	/
994	210112217	M	24	Oropharyngeal Swab	/	-	-	> 40	/
995	210112218	M	40	Oropharyngeal Swab	/	-	-	> 40	/
996	210112219	M	34	Oropharyngeal Swab	/	-	-	> 40	/
997	210112220	F	39	Oropharyngeal Swab	/	-	-	> 40	/
998	210112221	F	28	Oropharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thuringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

999	210112222	F	36	Oropharyngeal Swab	/	-	-	> 40	/
1000	210112223	M	41	Oropharyngeal Swab	/	-	-	> 40	/
1001	210112224	F	26	Oropharyngeal Swab	/	-	-	> 40	/
1002	210112225	M	39	Oropharyngeal Swab	/	-	-	> 40	/
1003	210112226	F	31	Oropharyngeal Swab	/	-	-	> 40	/
1004	210112227	F	20	Oropharyngeal Swab	/	-	-	> 40	/
1005	210112228	M	26	Oropharyngeal Swab	/	-	-	> 40	/
1006	210112229	F	28	Oropharyngeal Swab	/	-	-	> 40	/
1007	210112230	M	30	Oropharyngeal Swab	/	-	-	> 40	/
1008	210112231	M	30	Oropharyngeal Swab	/	-	-	> 40	/
1009	210112232	F	32	Oropharyngeal Swab	/	-	-	> 40	/
1010	210112233	F	44	Oropharyngeal Swab	/	-	-	> 40	/
1011	210112234	F	32	Oropharyngeal Swab	/	-	-	> 40	/
1012	210112235	F	42	Oropharyngeal Swab	/	-	-	> 40	/
1013	210112236	M	40	Oropharyngeal Swab	/	+	-	> 40	/
1014	210112237	M	32	Oropharyngeal Swab	/	-	-	> 40	/
1015	210112238	F	31	Oropharyngeal Swab	/	-	-	> 40	/
1016	210112239	F	38	Oropharyngeal Swab	/	-	-	> 40	/
1017	210112240	M	46	Oropharyngeal Swab	/	-	-	> 40	/
1018	210112241	M	43	Oropharyngeal Swab	/	-	-	> 40	/
1019	210112242	M	42	Oropharyngeal Swab	/	-	-	> 40	/
1020	210112243	F	32	Oropharyngeal Swab	/	-	-	> 40	/



MP Biomedicals Germany GmbH
 Thueringer Straße 15
 37269 Eschwege GERMANY
 mpbio.com

1021	210112244	F	53	Oropharyngeal Swab	/	-	-	> 40	/
1022	210112245	M	34	Oropharyngeal Swab	/	-	-	> 40	/
1023	210112246	F	29	Oropharyngeal Swab	/	-	-	> 40	/
1024	210112247	M	43	Oropharyngeal Swab	/	-	-	> 40	/
1025	210112248	M	24	Oropharyngeal Swab	/	-	-	> 40	/
1026	210112249	F	31	Oropharyngeal Swab	/	-	-	> 40	/
1027	210112250	F	22	Oropharyngeal Swab	/	-	-	> 40	/